

Dermapen HOME™

Start your microneedling journey with confidence

From setting up your device to accelerating results, discover everything you need to know about microneedling safely at home with our easy-to-follow 8-step guide.





Why Dermapen HOME™?

If you've already embarked on your microneedling journey in a clinic, it's time to get your hands on Dermapen HOME™.

With a little practice, you'll confidently perform 'top-up' treatments between clinic appointments—all from the comfort of your home. And, the results will be worth it!

When following the steps in this guide, Dermapen HOME™ is a safe and easy way to achieve skin that looks and feels amazing.

Endorsed by leading skin specialists, Dermapen HOME™ is shown to complement professional treatments for a wide range of skin concerns. From ageing skin to pigmentation, rosacea, acne, scars and more, it's one device to combat multiple skin conditions.



For step-by-step instructions, scan the QR code and watch the video. It makes using your Dermapen HOME™ easy.



01


Set Up Your Dermapen HOME™

It's easy to get started, but it's important to follow these quick steps before you begin:

- Install the DermapenWorld™ app.
- Create an account.
- Register your device to the account.

Your device comes with two DermapenWorld™ turbo-cut lithium batteries that are precharged, allowing you to perform treatments immediately.

TIP: Keep one battery on charge at all times as your backup (it takes 2 hours to fully charge).

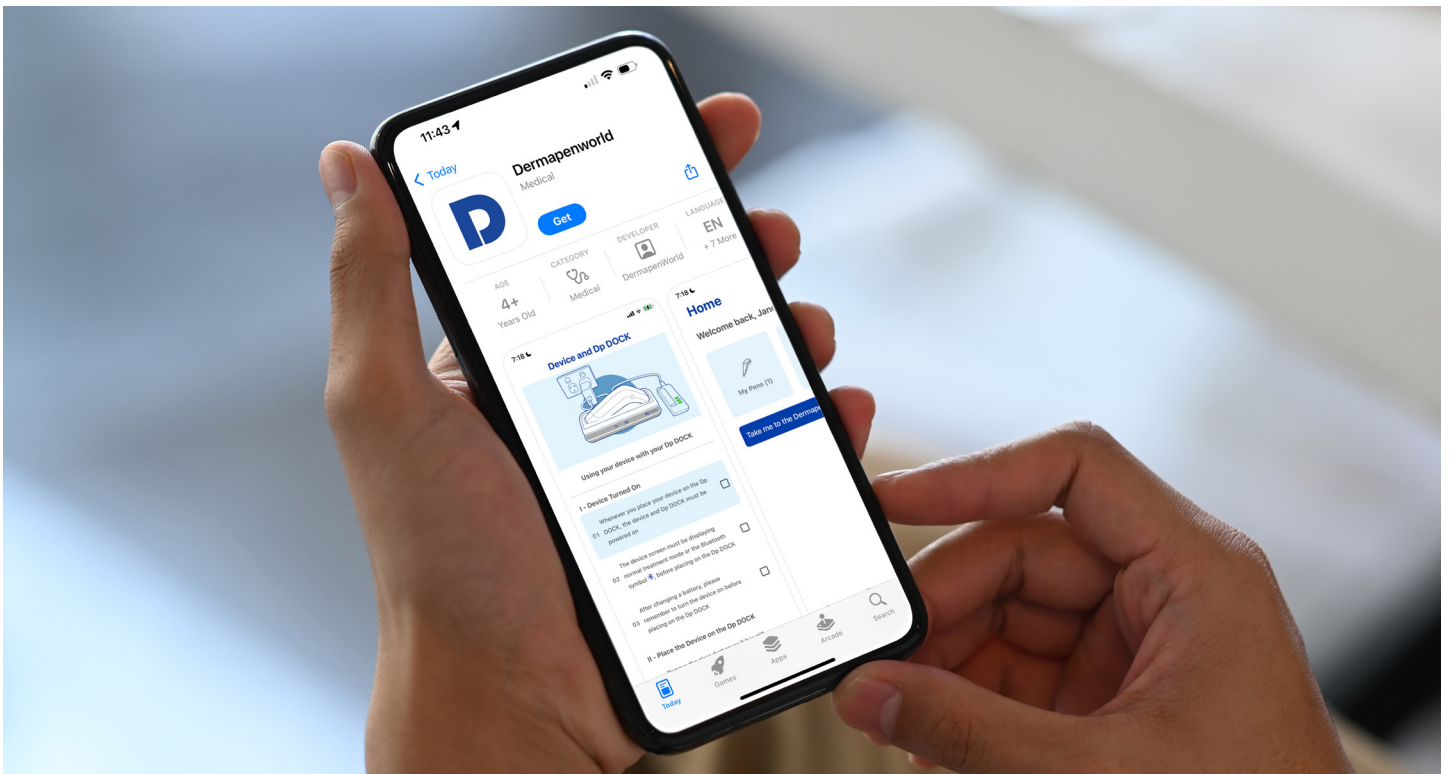

02

Choose Your Meso-Glide™ Serum



Microneedling is best performed with a serum or cream called a Meso-Glide™.

Dermapen HOME™ creates micro-channels, which deliver active ingredients of the Meso-Glide™ into the skin—up to 80% deeper than regular topical applications.

Applying a Meso-Glide™ will also help your device to glide over the skin for a comfortable experience.

Your Dermapen HOME™ device comes with the right amount of HYLA ACTIVE™ so you can start your first procedure immediately (2-3mls is all you need).

HYLA ACTIVE™ delivers deep hydration, comfort, and repair when used on its own, or you can mix it with the following serums.



CLR LOTION™
mixed 1:1 with HYLA ACTIVE™
Target acne and breakouts with proven ingredients while hydrating and balancing oil production.



BRITE LITE™
mixed 1:1 with HYLA ACTIVE™
Help fade pigmentation and dark spots and correct uneven skin tone—even dark undereye circles.



ANTIOXIDANT COCKTAIL™
mixed 1:1 with HYLA ACTIVE™
Combat signs of ageing and sun damage, including fine lines and wrinkles, thinning or sagging skin and a dull complexion.

03

Thoroughly Wash Hands and Face

Cleanliness is key to effective and safe home treatments.

Make sure your surface area and pen are clean and disinfected. Have fresh hand towels or gauze ready. Dispense and pre-mix your Meso-Glide™ in a clean bowl so that it's easy to apply during treatments.

We recommend Dp Dermaceuticals™ TRI-PHASE CLEANSER™. Ensure skin is dry before spritzing with CLINIPREP™ for clinical standard hygiene.

Depending on skin type, you may also use CLR FOAM CLEANSER™ or MICRO DERM EXFOLIANT™.

TIP: Cleanse hands with Dp Dermaceuticals™ HAND SANITISER™.



04

Attach Sterile Needle Cartridge

With clean hands, remove Needle Cartridge from sealed blister pack. Gently insert it into Dermapen HOME™.

Turn the cartridge a quarter-turn clockwise to lock in place. It should slot in easily—don't force, squeeze or overturn as this may damage the device or cartridge.

TIP: For your safety, needle cartridges are designed for single use. After the treatment, remove cartridge by gently rotating it in an anti-clockwise direction.

Using force to remove the cartridge may cause damage to the device.



05

Set Needle Level and Speed


Switch Dermapen HOME™ on by holding the Primary Button for **3 seconds**.

Remove the protective cap from the Needle Cartridge and discard.

Use the Level Buttons on the top of the device to increase (+) or decrease (-) the needle level.

Use the Speed Buttons on the side of the device to increase (+) or decrease (-) the needle speed.

Primary Button

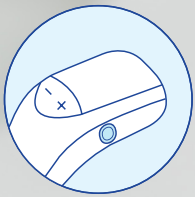


EXPERT TIP

High speed encourages faster results, and a higher needle setting may achieve visible results more quickly. However, if it isn't comfortable, don't worry. Start on the lowest level and dial-up to the maximum level according to your tolerance.

You can still achieve outstanding results by performing Dermapen HOME™ procedures at the minimum level every 7-10 days.


Speed Buttons

Level Buttons

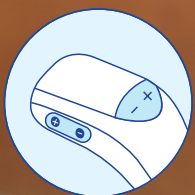

06

Apply Meso-Glide™ in Sections

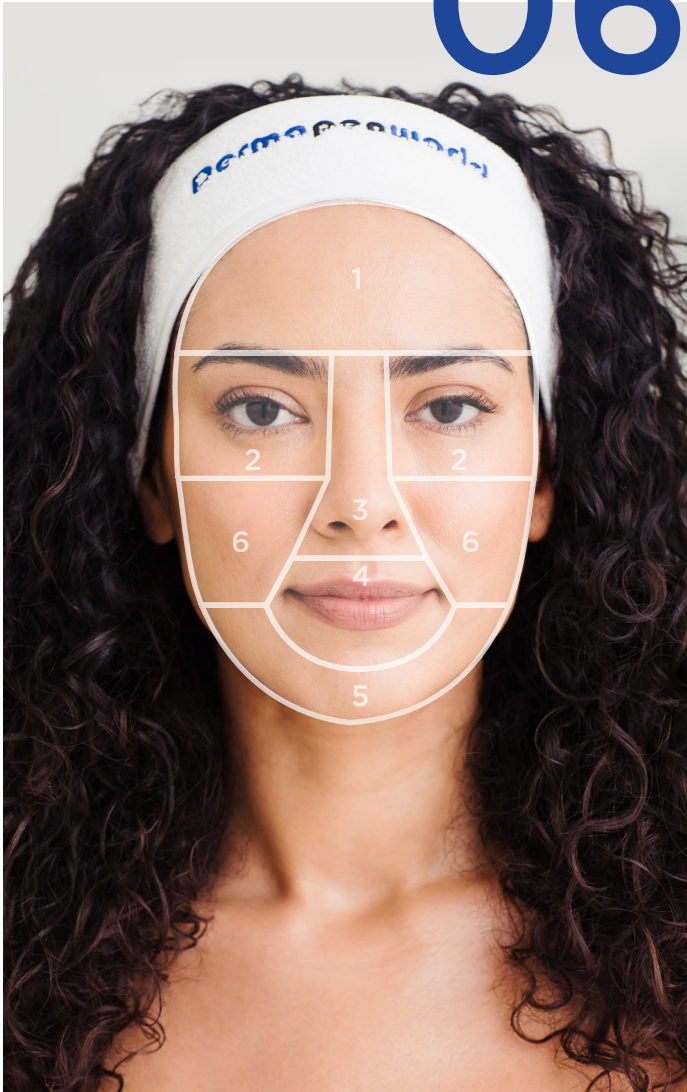

Start with one section of the face at a time and ensure it's completely covered with serum. As you progress through the treatment, liberally apply more serum to the next section. Your skin will soak up the nutrients and vitamins faster than ever!

Split the face into the following six areas to ensure a thorough treatment. Some areas may require the needle level or speed adjusted to your comfort level.

1. **Forehead:** May feel slightly sensitive, but imagine those frown lines disappearing.
2. **Eye area:** For safety, avoid needling too close to the eyes. Areas outside the orbital area, such as crow's feet or dark circles below this delicate area, may be treated with caution.
3. **Nose:** Another sensitive area that's often impacted by enlarged pores, blackheads, and sun damage.
4. **Upper lip and lips:** Yes, you can plump your lips, too! Lines can easily form on the upper lip when it's overlooked during skincare routines.
5. **Chin/jawline:** Minimise pores and breakouts in this area and tighten jowls.
6. **Cheeks:** Often where signs of sagging or loose skin, uneven skin tone or texture, you can go deeper and faster on this area.

Microneedling beyond your face:

- **Neck:** Target lines and drooping skin.
- **Décolleté:** An often-neglected area, which shows age through thinning and sun damage.
- **Ears:** Behind the ears is often a tell-tale sign of age, and the ear lobes also begin to "sag" as we age (we call this the "Chandlear Treatment").
- **Hands:** Reduce pigmentation and plump thinning or dehydrated skin.
- **Knees:** Give your legs a make-over by rejuvenating sagging knees.
- **Feet:** Keep your feet looking their best.

And anywhere else you feel needs a boost of rejuvenating collagen.

07

Start Microneedling


Glide the scalloped rim of the Needle Cartridge using an even pressure over the skin to ensure constant needle contact and distribution of the Meso-Glide™. Use stripe-like motions, lifting the needles between each stripe to avoid unnecessary reinsertions.

Complete 4 motions


1. Up and down




2. Side to side

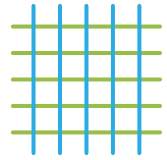

3. Diagonal direction

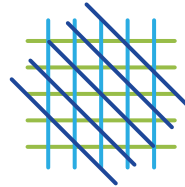


4. Opposite diagonal direction

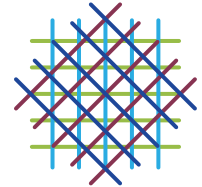

Think of the Union Jack flag as your guide.

Repeat within the treatment area 3-5 times before moving onto the next zone. You can adjust speed and level as you go.

Apply more Meso-Glide™ throughout your treatment.

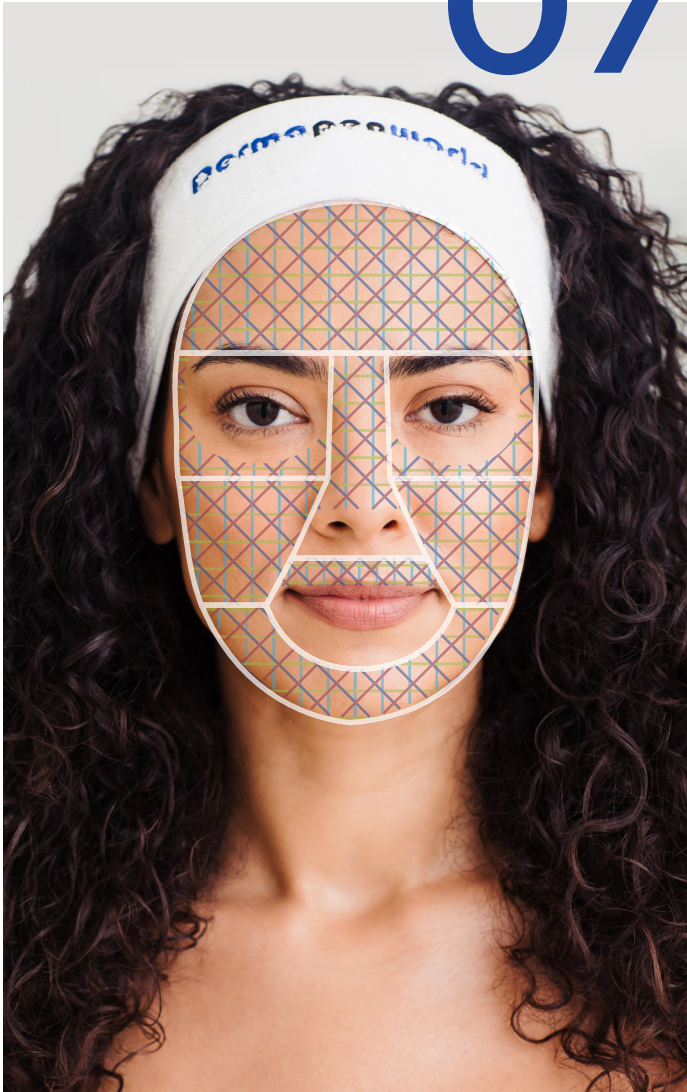

08

Accelerate Healing

If you experience inflammation and sensitivity following your treatment, you're not alone. We recommend you calm and soothe your skin with a Dp Dermaceuticals™ Hyla Active™ 3D Sculptured Mask. Apply for 10-30 minutes.

Follow with Dp Dermaceuticals™ Vitamin Rich Repair™.

As a final step, apply Dp Dermaceuticals™ Cover Recover™ to camouflage any redness, protect from UV rays, and further soothe the skin.



Maximise Results

By incorporating LED light therapy into your microneedling program you can ramp up your results.

LED light therapy is proven to:

- reduce inflammation
- assist skin to heal
- increase the absorption of skincare products.

Use a Dp Dermaceuticals L.E.D. red and near-infrared light mask immediately after microneedling and in between treatments for beautiful glowing skin.

For even more benefits, pair L.E.D. FACE with LUMAFUSE™ HYDROGEL SHEET MASK. The light helps the actives reach deeper into the skin, and the jelly-like material of the mask helps the light to penetrate deeper too. Enjoy smooth, radiant skin even faster.



Dermapenworld™

www.DermapenWorld.com