

# Dermapenworld™

OUR SYNERGISTIC SKINCARE SOLUTIONS

Dp Dermaceuticals L.E.D.:  
Getting Visible Results At Home



Home protocols are an important part of our synergised solutions to achieve and keep the results you want.

Follow this guide to perform your at-home L.E.D. Light Therapy treatments with complete confidence.



## 01 Enjoy Clinically Proven Technology

Almost 100% of people surveyed said they saw visible improvements when following the recommended treatment programme with Dp Dermaceuticals L.E.D.

Set aside 10 minutes, 3-4 times a week, for the next 4 weeks (and make sure to take your before pic now!).

## 02 Charge L.E.D.

### L.E.D. FACE AND L.E.D. NECK & DÉCOLLETÉ:

Connect the power adaptor to the power cable, and plug into a USB-C charging port. This could be a wall socket, laptop or any power outlet compatible with USB-C.

Check the three dots on the power adaptor light up one after another. This shows L.E.D. is charging successfully.

### L.E.D. HAND:

Connect L.E.D. hand directly to the power cable and plug into a USB-C charging port.





## 03 Unplug the Power Cable

A safety feature of L.E.D. is that it won't turn on while connected to the mains power supply - this means you can safely wear your L.E.D. anywhere, even in the bath! When all three dots are illuminated continuously, your L.E.D. is fully charged. Disconnect the power cable.

## 04 Adjust Straps

Dp Dermaceuticals L.E.D.s are fully adjustable and designed to stay in place while sitting, lying or moving around. Try the straps at different lengths and positions until you get the right fit for your face, head, neck, décolleté and wrists.

### **EXPERT TIP:**

L.E.D. FACE straps are best positioned higher on the crown of the head to ensure gravity doesn't interrupt your treatment with any slipping. The device should be comfortable and secure no matter what else you're doing.





## 05 Apply LumaFuse Hydrogel Mask or Serum

Intensify L.E.D. treatments with the LumaFuse Hydrogel Mask, specifically designed to deliver hydrators, collagen, antioxidants, plant extracts and vitamins in combination with the maximum penetration of red and near infrared light rays absorbed into the skin.

Alternatively, smooth Hyla Active or one of our correcting serums to skin under the L.E.D. and allow to absorb naturally during the treatment.



### HYLA ACTIVE

Potent breakthrough Hyaluronic Acid (HA) formulation delivers biologically active HA, Zinc and Vitamin B5 deeper into skin, providing deep hydration, comfort and repair.



### ANTIOXIDANT COCKTAIL

Combat signs of Ageing and Sun Damage including fine lines and wrinkles, thinning or sagging skin and brighten dull complexions.




### BRITE LITE

Dual-action wonder serum to fade pigmentation and dark spots, and fix uneven skin tone – even dark undereye circles.



### CLR LOTION

Treat Acne & Breakouts with proven ingredients. Remove impurities, hydrate and balance oil production.

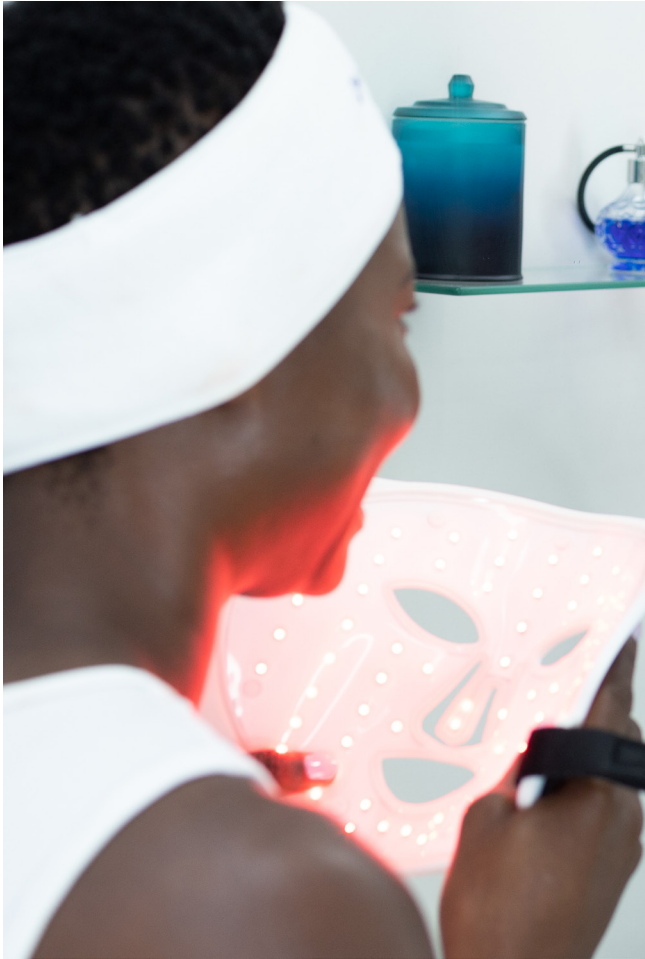

## 06 Start Treatment

When the L.E.D. is fitted comfortably, press the ON button on the power adaptor for 3 seconds or until the L.E.D. lights up.

Enjoy the soothing effect as it boosts your serotonin levels, as well as treating the skin!

The L.E.D. automatically turns off after 10 minutes but double dosing is completely safe if you're not ready to end the treatment.

Please note, the light may be very bright initially when treating the face. It is proven completely safe and will not damage the eyes in any way. You can wear a sleeping mask underneath if the light is uncomfortable.

## 07 Treatment Programme

To achieve the clinically proven results and experience plumped skin with visibly diminished fine lines, reduced pigmentation and smooth, glowing skin, repeat the treatment 3-4 times every week for 4 weeks.

Treat yourself before bedtime and enjoy increased serotonin levels and a deeper sleep.





#### EXPERT TIP:

L.E.D. Light Therapy has been shown to decrease pain too. Use regularly to relieve aches and inflammation from injuries, surgery or treatments such as microneedling procedures. Downtime following Dermapen Treatments is significantly reduced and results accelerated by applying L.E.D. immediately after.

Dermapen Treatments are available in clinics or you can safely perform treatments at home yourself with Dermapen HOME.

## 08 Clean L.E.D.

Always clean your L.E.D. with CLINIPREP before and after treatments. Stubborn stains can be removed with 70% alcohol.

Store L.E.D. flat in the dust bag provided.

