

Liquid Scope HelpOR™ II


Our vision is your clear vision.



Simple, instant and effective.

The Symmetry Surgical® **Liquid Scope HelpOR™ II** is a cost-effective and comprehensive scope management system that enhances a surgeon's vision during laparoscopic and robotic surgeries. The scope management system features an extended three-button Liquid Scope WarmOR to de-fog the scope lens as well as a lint-free microfiber cloth and trocar swabs to efficiently clean the scope lens and remove debris from the cannula.

Scopes are defogged and cleansed quickly and efficiently without costly capital equipment or electrical connections.

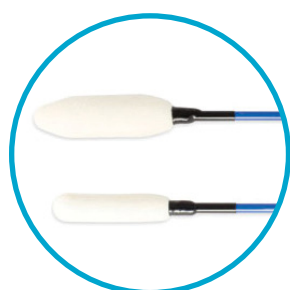

Trocar Cleaning Swabs



Liquid Scope HelpOR




Long Scope WarmOR



A SIMPLE AND INSTANT HEAT SOURCE


- Activates quickly and intuitively
- Reaches peak temperature (122°F | 50°C) in 5 minutes
- Retains optimal temperature for up to one hour
- Improves OR safety by eliminating fire, burn, and hot water hazards
- Accommodates all sizes and lengths of endoscopes
- Does not require batteries or electrical connections
- Is an affordable, disposable, and sterile device
- Is latex, DEHP, and alcohol-free.

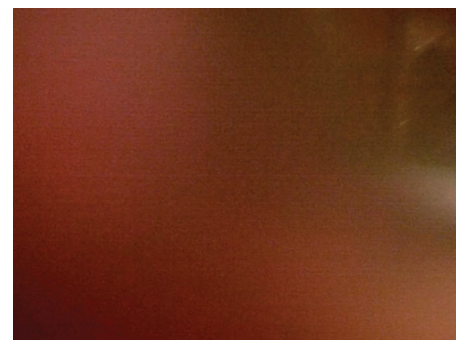

Fogged view when not using Liquid Scope HelpOR

CLEAN SCOPES AND TROCARS

- Includes a non-abrasive and soft microfiber cleaning cloth (15 x 15cm) that is gentle on surfaces without streaking, marring, or scattering
- Includes two trocar cleaning swabs (for 5-8mm and 10-12mm trocars) with radiopaque flexible shafts for cleaning in tight spaces
- Includes trocar cleaning swabs with soft, atraumatic foam tips that remove the blood and other fluids from the trocar and cannula without damaging the valves.

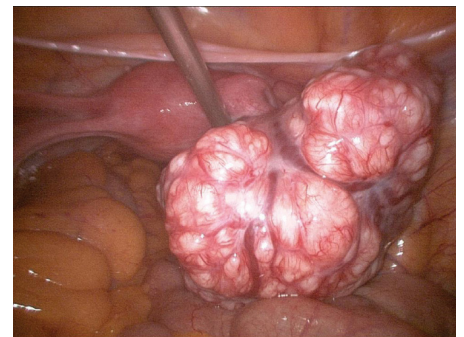

Clear view when using Liquid Scope HelpOR

REF / Order #	Description	UOM
LSH-02	Liquid Scope WarmOR™ long (three buttons), two trocar swabs (5-8mm and 10-12mm), one microfiber cloth	Box of 10, sterile

Contact a Surgical Specialist to learn more about the Minimally Invasive Surgery portfolio.


symmetry surgical® | **Minimally Invasive Surgery**
 customerservice@symmetrysurgical.com • 800.251.3000 • www.symmetrysurgical.com

08.19.2022 v1