

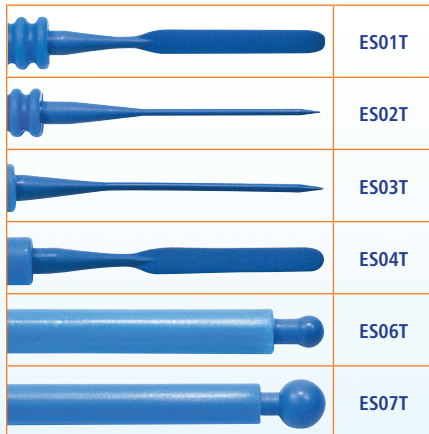


RESISTICK II CAN SAVE TIME & MONEY! RESISTS ESCHAR BUILD-UP & CAN SPEED UP PROCEDURES

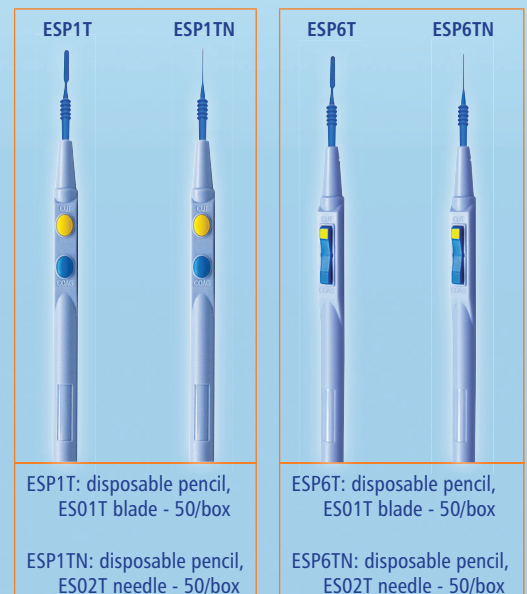

Resistick II™ coated blade, ball and needle electrodes are expertly coated with PTFE (polytetrafluoroethylene). They move more effectively through tissue because they resist eschar build-up making them much more efficient during surgical procedures. Unlike the common stainless steel electrodes with a propensity to accumulate the discarded blood and tissue, Resistick II™ coated electrodes simply wipe clean. Just use a damp sponge or cloth to almost effortlessly clean the debris off the Resistick II™ electrodes. Standard stainless steel electrodes are unquestionably much more difficult to clean. As they say, "Time is Money"!

To guarantee the consistency and thickness of the Resistick II™ coating, Bovie® utilizes an Ionic Backscatter measuring system which helps promote accurate and consistent coating thickness. By controlling the thickness of the coating, Resistick II™ by Bovie® gives the surgeon the optimal performance in each of the 20 different coated electrode styles.

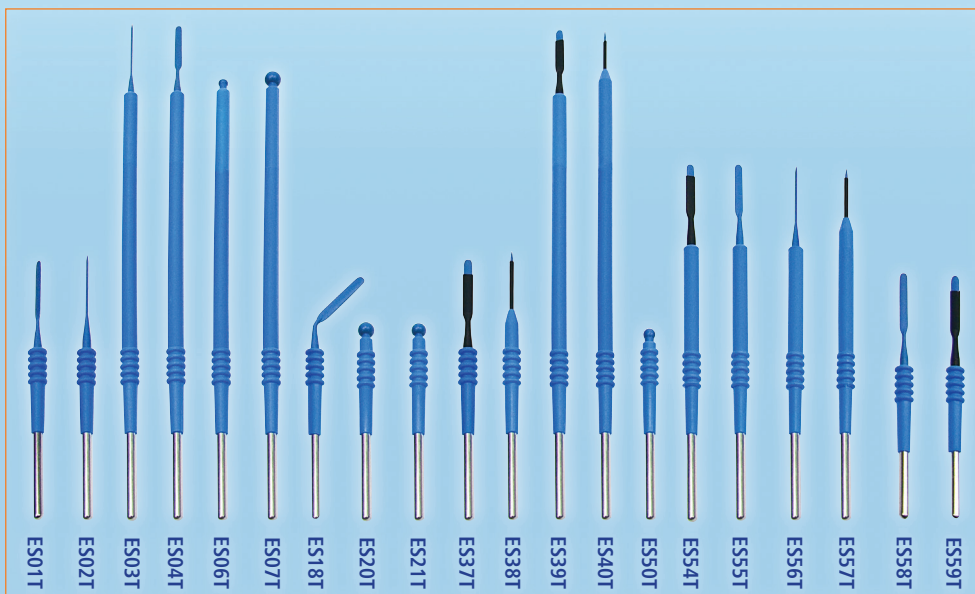

BOVIE®	DESCRIPTION	VALLEYLAB™ (EDGE™)	CONMED® (ULTRA CLEAN®)	MEGADYNE® (E-Z CLEAN®)
ES01T	RESISTICK II COATED STANDARD BLADE ELECTRODE 2.75" (6.98 cm)	E1425X	N/A	0012A
ES02T	RESISTICK II COATED STANDARD NEEDLE ELECTRODE 2.75" (6.98 cm)	E1452	139102	0013
ES03T	RESISTICK II COATED EXTENDED NEEDLE ELECTRODE 6" (15.24 cm)	E1452-6	N/A	0016
ES04T	RESISTICK II COATED EXTENDED BLADE ELECTRODE 6" (15.24 cm)	E1450-6	139107	N/A
ES06T	RESISTICK II COATED EXTENDED 3MM BALL ELECTRODE 5" (12.70 cm)	N/A	N/A	N/A
ES07T	RESISTICK II COATED EXTENDED 5MM BALL ELECTRODE 5" (12.70 cm)	N/A	N/A	0009
ES18T	RESISTICK II COATED ANGLED BLADE 45° ELECTRODE 2.75" (6.98 cm)	N/A	N/A	N/A
ES20T	RESISTICK II COATED STANDARD 5mm BALL ELECTRODE 2" (5.08 cm)	N/A	N/A	0015
ES21T	RESISTICK II COATED STANDARD 4mm BALL ELECTRODE 2" (5.08 cm)	N/A	N/A	N/A
ES37T	RESISTICK II COATED MODIFIED BLADE ELECTRODE 2.75" (6.98 cm)	E1455	139104EXT	0012AM
ES38T	RESISTICK II COATED MODIFIED NEEDLE ELECTRODE 2.75" (6.98 cm)	E1465	139105EXT	0013M
ES39T	RESISTICK II COATED EXTENDED MODIFIED BLADE ELECTRODE 6" (15.24 cm)	E1455-6	139110EXT	N/A
ES40T	RESISTICK II COATED EXTENDED MODIFIED NEEDLE ELECTRODE 6" (15.24 cm)	N/A	139025EXT	0016M
ES50T	RESISTICK II COATED STANDARD 3mm BALL ELECTRODE 2" (5.08 cm)	N/A	N/A	N/A
ES54T	RESISTICK II COATED EXTENDED MODIFIED BLADE ELECTRODE 4" (10.16 cm)	E1455-4	139112EXT	0014AM
ES55T	RESISTICK II COATED EXTENDED BLADE ELECTRODE 4" (10.16 cm)	E1450-4	139112	0014A
ES56T	RESISTICK II COATED EXTENDED NEEDLE ELECTRODE 4" (10.16 cm)	N/A	N/A	0016A
ES57T	RESISTICK II COATED EXTENDED MODIFIED NEEDLE ELECTRODE 4" (10.16 cm)	E1465-4	N/A	0016AM
ES58T	RESISTICK II COATED BLADE ELECTRODE 2.5" (6.35 cm)	E1450G	139100	0012
ES59T	RESISTICK II COATED MODIFIED BLADE ELECTRODE 2.5" (6.35 cm)	N/A	N/A	0012M

Bovie is a registered trademark of Bovie Medical Corporation.
 Resistick is a trademark of Bovie Medical Corporation.
 Valleylab and Edge are trademarks of Covidien.
 UltraClean is a registered trademark of ConMed Electrosurgery.
 E-Z Clean is a registered trademark of Megadyne.

