





### WILLIAM T. BOVIE - THE INVENTOR OF ELECTROSURGERY

The first use of the Bovie® electrosurgical generator in an operating room was on October 1, 1926, at Peter Bent Brigham Hospital in Boston, Massachusetts. A few days earlier (September 28, 1926) Dr. Harvey Cushing had been unsuccessful in removing a mass from a patient's head due to its highly vascular nature. In Dr. Cushing's notes from October 1, "...with Dr. Bovie's help I proceeded to take off most satisfactorily the remaining portion of tumor with practically none of the bleeding which was occasioned in the preceding operation."

Almost a century later, Symmetry Surgical® continues to improve the work of Dr. Bovie. This catalog represents our Energy Solution portfolio, including Bovie electrosurgical generators, high-frequency desiccators, monopolar and Olsen® bipolar instrumentation, smoke evacuation, lighting, cauteries, colposcopes and more. The work continues today, but our sincere thanks go out to Dr. Cushing and Dr. Bovie for their incredible invention.

The Bovie® electrosurgery portfolio is one of Symmetry Surgical's valued brands. Symmetry Surgical is a leading healthcare company that develops, manufactures, and markets high-quality surgical instrumentation, energy and smoke evacuation products that meet clinician needs to positively impact patients and the economics of surgical care. Our broad portfolio of single-use, limited-use and reusable surgical instruments are used in every surgical specialty in sites of care around the world.

### **Symmetry Surgical**

Antioch, TN 37013, USA
U.S. Phone 1-800-251-3000
symmetrysurgical.com
boviemedical.com
customerservice@symmetrysurgical.com





### **CONTENTS**

ELECTRUSURGERY	6 - 29
<ul> <li>OR   PRO 300</li> <li>Surgi-Center   PRO</li> <li>Specialist   PRO</li> <li>Specialist   PRO-G</li> <li>Bantam   PRO</li> <li>Bantam   PRO-G</li> <li>DERM 942</li> <li>DERM 102</li> <li>Electrosurgery Pencils</li> <li>Suction Coagulators</li> <li>Disposable Electrodes, Blades, Balls, Needles, Dermal Tips</li> <li>Resistick II<sup>™</sup> Coated Electrodes</li> <li>Olsen<sup>®</sup> and Bovie<sup>®</sup> Microdissection Needles</li> <li>Laparoscopic Electrodes</li> <li>Disposable Loops</li> <li>Tungsten SuperCut<sup>™</sup> Needles &amp; Loop, Arthroscopic Electrodes</li> <li>Reusable Electrodes, Blades, Balls, Needles</li> <li>Reusable Bipolar Forceps</li> <li>Olsen<sup>®</sup> Bipolar Forceps</li> <li>Mobile Stands, Footswitches, Wall Mount Kits</li> <li>Return Electrodes, Active Cable, Cords, Adaptors</li> </ul>	7 8 9 10 11 12 13 14 15 16 17 18 19 20 21 22 23 24 25 26 - 27 28 29
SMOKE EVACUATION	30 - 34
<ul> <li>Bovie Smoke Shark II</li> <li>Bovie Smoke Shark III</li> <li>Smoke Evacuation Pencils and Adaptors</li> <li>Smoke Evacuation Accessories</li> </ul>	31 32 33 34
MEDICAL LIGHTS	35 - 39
<ul> <li>MI 150 and MI 300 LED Lights</li> <li>MI 550 LED Lights</li> <li>MI 750 LED Lights</li> <li>MI 1000 LED Lights</li> <li>System Two LED Lights</li> </ul>	35 36 37 38 39

COLOSCOPES AND ACUDOP™ II	40 - 44
<ul> <li>Colpo-Master™ I</li> <li>Colpo-Master™ II</li> <li>Olsen® Vaginal Speculum with Smoke Port</li> <li>AcuDop™ II</li> </ul>	41 42 43 44
CAUTERIES	45 - 49
<ul> <li>High-Temperature Cauteries</li> <li>Low-Temperature Cauteries</li> <li>Change-A-Tip™ Cauteries</li> <li>Replacement Tips</li> </ul>	46 47 48 49
SPECIALTY LIGHTS	50 - 52
<ul><li>Flexible Lights, Stylets</li><li>Specialty Lights, Penlights, UV Woods</li></ul>	51 52
OTHER PRODUCTS	53 - 54
<ul> <li>DermaScout<sup>™</sup> II</li> <li>Burr Handles, Burr Tips, Eye Bubble</li> <li>Neuro-Pulse<sup>™</sup></li> </ul>	53 54 54

### **ELECTROSURGERY**



### **BOVIE® OR | PRO 300 ELECTROSURGICAL GENERATOR**

The Bovie® OR | PRO 300 system was designed in conjunction with O.R. clinicians to provide a safe, user-friendly generator to meet all operating room procedure performance demands. This innovative energy-based device offers 13 modes of surgical energy delivery. Special modalities include the advanced Bovie® Bipolar, Auto Bipolar options, as well as special Laparoscopic modes that enable surgeons to select the right generator to meet their requirements for a variety of surgical applications and preferences.



## RECOMMENDED ACCESSORIES † BV-1253B BV-1254B ESRE A821 A1255A A827V ESP1, ESP6 or ESPR2 BV-IDS-CS2

### Style # A3350

### **▼** Products Included

**A3350** (1) electrosurgical generator 110 VAC Hospital-grade power cord\* 220 VAC Hospital-grade power cord\* User's guide (MC-55-230-001)\*

### **Dimensions**

Height: 16.5cm (6.5") Width: 37.5cm (14.75") Depth: 46cm (18.1") Weight: <9.07kg (<20lbs)

### **Output Power**

Cut I & II: 300 watts
Blend (1-4): 200 watts
Pinpoint Coagulation: 120 watts
Spray Coagulation: 120 watts
Gentle Coagulation: 120 watts
Macro Bipolar: 150 watts
Micro Bipolar: 80 watts
Standard Bipolar: 50 watts
Bovie\* Bipolar: 225
Auto Bipolar: 80

### **Output Frequency**

490 kHz ± 5 kHz

### Line Voltage

100 - 240 VAC

### **Line Frequency**

50 - 60 Hz

- \* Not Pictured
- † Sold Separately

### **BOVIE® SURGI-CENTER | PRO ELECTROSURGICAL GENERATOR**

The Bovie® Surgi-Center | PRO Electrosurgical Generator is commonly used for cutting, coagulation and ablation of tissue in general, gynecologic, orthopedic, ENT and urological procedures performed in an operating suite or procedure room. The unit was designed to provide an electrosurgery solution capable of most operating room demands in an affordable 200 watt generator. The unit offers 12 options of energy delivery for a variety of surgical application; two cut, four blend, three coagulation and three bipolar modalities.



# ★ RECOMMENDED ACCESSORIES † BV-1253B BV-1254B ESRE A821 A1255A A827V ESP1, ESP6 or ESPR2 BV-IDS-CS2

### Style # A2350

### **▼ Products Included**

**A2350** (1) electrosurgical generator 110 VAC Hospital-grade power cord\* 220 VAC Hospital-grade power cord\* User's guide (**MC-55-232-001**)\*

### **Dimensions**

Height: 16.5cm (6.5") Width: 37.5cm (14.75") Depth: 46cm (18.1") Weight: <9.07kg (<20 lbs)

### **Output Power**

Cut I & II: 200 watts Blend (1-4): 200 watts Pinpoint Coagulation: 120 watts Spray Coagulation: 120 watts Gentle Coagulation: 120 watts Macro Bipolar: 80 watts Micro Bipolar: 80 watts Standard Bipolar: 50 watts

### **Output Frequency**

490 kHz ± 5 kHz

### **Line Voltage**

100 - 240 VAC

### **Line Frequency**

50 - 60 Hz

- \* Not Pictured
- <sup>†</sup> Sold Separately

### **SPECIALIST | PRO HIGH-FREQUENCY ELECTROSURGICAL GENERATOR**

The Specialist | PRO is the affordable solution for the busy physician office or specialty facility. Capable of most monopolar and bipolar surgical procedures and packed with dependable safety features, the Specialist | PRO satisfies the physician's surgical demands with safety, reliability and convenience. When paired with the Bovie® Smoke Shark® II Smoke Evacuator, it is a total system solution for procedures and facilities that require smoke evacuation.



### **▼ PRODUCTS INCLUDED** ES01, ES02, ES20 **FSRF** A1252C

### **★ RECOMMENDED ACCESSORIES** <sup>†</sup> BV-1253B A1255A

### **Style # A1250S**

### **▼** Products Included

A1250S (1) electrosurgical generator ESP1 (1) sterile disposable pencil (1) blade sterile disposable electrode **ES01** ES02 (1) needle sterile disposable electrode **ES20** (1) ball sterile disposable electrode **ESRE** (5) disposable split return electrodes A1252C (2) reusable grounding cords 110 VAC Hospital-grade power cord\* 220 VAC Hospital-grade power cord\*

### **Dimensions**

Height: 15.2cm (6") Width: 26cm (10.25") Depth: 30.5cm (12") Weight: < 6.5kg (< 14 lbs.)

### **Output Power**

Cut: 120 watts Blend: 90 watts Coagulation: 80 watts Fulguration: 40 watts Bipolar: 30 watts

### **Output Frequency**

350 – 600 kHz

### Line Voltage

100-240 VAC

### **Line Frequency**

50 - 60 Hz

- \* Not Pictured
- <sup>†</sup> Sold Separately

### SPECIALIST | PRO-G OB/GYN TOTAL SYSTEM SOLUTION

The A1250S-G combines the Specialist | PRO package with a Bovie® Smoke Shark® II smoke evacuation system and the latest mobile stand technology to produce a compact, mobile system for the office or specialty facility. As an added benefit, the A1250S-G includes a loop electrode kit (five of our most popular loops), five extended ball electrodes and a smoke evacuator hose reducer fitting for attaching to a vaginal speculum.



### **Style # A1250S-G**

### **▼** Products Included

### **Electrosurgery equipment**

**A1250S** (1) electrosurgical generator

**ESP1** (1) disposable pencil

ES01 (1) blade sterile disposable electrode ES02 (1) needle sterile disposable electrode

ES20 (1) ball sterile disposable electrode
ESRE (5) disposable split return electrodes

**A1252C** (2) reusable grounding cords **ES07** (5) extended ball electrodes 110 VAC Hospital-grade power cord\*

220 VAC Hospital-grade power cord\* User's guide (MC-55-237-001)\*

### **ESLK** electrosurgical loop kit

ES09 (1) 10mm x 10mm loop ES11 (1) 15mm x 10mm loop ES12 (1) 20mm x 8mm loop ES13 (1) 20mm x 15mm loop ES16 (1) 10mm x 8mm square

### **Mobile stand**

**ESMS2** (1) mobile stand for A1250S

### Smoke evacuation system

**SE02** (1) Smoke Shark® II with pneumatic

footswitch

**SF35** (1) 35-hour filter

**786T** (1) 7/8" (22.2mm) tube x 6' (182.9cm)

ong

**SERF** (1) reducer fitting (attaches to vaginal

speculum)

110 VAC Hospital-grade power cord\*

User's guide (MC-23088)\*

### **BOVIE® BANTAM | PRO**

The Bovie® Bantam | PRO is a high frequency desiccator that also has cutting ability. The unit provides 50 watts of cutting power and coagulation modes — plenty of power for most procedures performed in the doctor's office. Power is adjustable in 1/10th watt increments from 0 to 10 watts, for the most accurate output available.



# ▼ PRODUCTS INCLUDED A902 A804, A805, A806, A807, A806DE, A807DE ES01, ES02, ES20 A837 A910 A910ST

### **Style # A952**

### **▼** Products Included

**A952** (1) electrosurgical generator

**A902** (1) handpiece

**A804** (5) non-sterile sharp dermal tips **A805** (2) sterile sharp dermal tips

**A806** (5) non-sterile blunt dermal tips **A807** (2) sterile blunt dermal tips

**A806DE**(2) non-sterile DERM-Elite™

Premium blunt dermal tips

**A807DE**(2) sterile DERM-Elite Premium blunt dermal tips

**ES01** (3) blade sterile disposable electrodes

**ES02** (1) needle sterile disposable electrodes

**ES20** (1) ball sterile disposable electrodes

**A910** (2) non-sterile disposable handpiece sheaths

**A910ST**(2) sterile disposable handpiece sheaths

A1252C (2) reusable grounding cord\*

**ESRE** (5) disposable split grounding pads\*

**A837** (1) wall mount kit

110 VAC Hospital-grade power cord\* User's guide (MC-55-238-001)\*

### **Dimensions**

Height: 188 mm (7.40 in.) Width: 228 mm (8.98 in.) Depth: 105 mm (4.13 in.) Weight: < 2.26 kg (< 5 lbs.)

### **Output Power**

Cut: 50 watts Blend: 50 watts Coagulation: 50 watts Fulguration: 40 watts Bipolar: 40 watts Micro Bipolar: 40 watts

Output Frequency: 304-451 kHz

Line Frequency: 50-60 Hz

Line Voltage: 100-240 VAC

### **BOVIE® BANTAM | PRO A952-G**

The A952-G combines the Bovie® Bantam | PRO A952 with the Bovie® Smoke Shark® II smoke evacuation system and a mobile stand to produce a compact, complete mobile system for the office. As an added benefit, the A952-G includes a loop electrode kit (five of our most popular loops), five extended ball electrodes and a smoke evacuation hose reducer that attaches to a vaginal speculum.



### **Style # A952-G**

### **▼** Products Included

Electrosur	gery	equipment
4050	14	\ I I

**A952** (1) electrosurgical generator

A902 (1) handpiece

A804 (5) non-sterile sharp dermal tips
A805 (2) sterile sharp dermal tips
A806 (5 non-sterile blunt dermal tips
A807 (2) sterile blunt dermal tips
(2) non-sterile DERM-Elite™
Premium blunt dermal tips

**A807DE** (2 sterile DERM-Elite Premium

blunt dermal tips

**ES01** (3) blade sterile disposable

electrodes

**ES02** (1) needle sterile disposable

electrodes

**ES20** (1 ball sterile disposable

electrodes

**A910** (2) non-sterile disposable

handpiece sheaths

**A910ST** (2) sterile disposable handpiece

sheaths

A1252C (2) reusable grounding cords\* (5) disposable split grounding

pads\*

**A837** (1) wall mount kit

**ES07** (5) extended ball electrodes 110 VAC Hospital-grade power cord\* User's guide (MC-55-238-001)\*

### **ESLK** electrosurgical loop kit

ES09 (1) 10mm x 10mm loop ES11 (1) 15mm x 10mm loop ES12 (1) 20mm x 8mm loop ES13 (1) 20mm x 15mm loop ES16 (1) 10mm x 8mm square

### **A812-C Mobile stand**

A812-E (1) mobile stand for A952 A812-E (1) height extender for A812 A808-T (1) top tray and clamp for A812 A812-BT (1) bottom tray for A812

Smoke evacuation system
SE02 (1) Smoke Shark® II with
pneumatic footswitch

**SF35** (1) 35-hour filter

**786T** (1) 7/8" (22.2mm) tube x 6'

(182.9cm) long

**SERF** (1) reducer fitting (attaches to

vaginal speculum)

110 VAC Hospital-grade power cord\* User's guide (MC-23088)\*

### **BOVIE® DERM 942**

The DERM 942 High Frequency Desiccator and its patent pending handpiece are designed and engineered to be reliable and durable. The unit's life is prolonged by incorporating the latest in chip technology and reducing the number of internal parts, thereby generating less internal heat. Post-procedure and long-term clean-up are simplified by its non-absorbing, non-staining, high impact plastic case. Power is adjustable in 1/10th watt increments in the first 10 watts, in monopolar and bipolar. The Bovie® DERM 942 is a ground reference unit, which means no grounding pad is required.



### ▶ PRODUCTS INCLUDED A902 A804, A805, A806, A807, A806DE, A807DE A837 A910

### **Style # A942**

### **▼** Products Included

A942 (1) high frequency desiccator

A902 (1) handpiece

**A804** (5) non-sterile sharp dermal tips **A805** (2) sterile sharp dermal tips **A806** (5) non-sterile blunt dermal tips

**A807** (2) sterile blunt dermal tips

**A806DE** (2) non-sterile DERM-Elite<sup>™</sup>

Premium blunt dermal tips

**A807DE** (2) sterile DERM-Elite Premium

blunt dermal tips

**A910** (2) non-sterile disposable

handpiece sheaths

**A837** (1) wall mount kit,

110 VAC Hospital-grade power cord\*
User's guide (**MC-55-239-001**)\*

### **Dimensions**

Height: 188mm (7.40") Width: 228mm (8.98") Depth: 105mm (4.13") Weight: < 2.26kg (< 5lbs.)

### **Output Power**

To 40 watts (0-10 watts in 1/10th watt increments on both monopolar and bipolar)

### Line Voltage

100-240 VAC

Output Frequency: 328-492 kHz

Line Frequency: 50-60 Hz

Warranty: FOUR YEARS

\*Not Pictured

### **DERM 102 HIGH FREQUENCY DESSICATOR**

The DERM 102 is an economical, 10 watt high frequency desiccator, thoughtfully designed to give practitioners the ability to perform in-office, minor skin procedures such as the removal of skin tags, benign lesions and premalignant lesions. The DERM 102 High Frequency Desiccator offers ten watts of monopolar and bipolar power, adjustable in 1/10th watt increments, to desiccate and fulgurate tissue. The power can be adjusted using the 3-Button Power Control Handpiece (A902) or using the rotary dial on the unit. The unit features a large, bright, LED display for excellent visibility across the procedure room.



### Style # DERM 102

### **▼** Products Included

**DERM 102** (1) high frequency desiccator

A902 (1) handpiece

A910 (2) non-sterile, disposable handpiece sheaths
A804 (5) non-sterile single-use sharp dermal tips
A805 (2) sterile single-use sharp dermal tips
A806 (5) non-sterile single-use blunt dermal tips
A807 (2) sterile single-use blunt dermal tips

**A806DE** (2) non-sterile DERM-Elite™ Premium blunt dermal tips **A807DE** (2) sterile DERM-Elite Premium blunt dermal tips

**A837** (1) wall mount kit

110 VAC Hospital-grade power cord\* User's quide CD (**MC-55-229-001**)\*

### **Dimensions**

Height: 208mm (8.2") Width: 184mm (7.2") Depth: 69mm (2.7") Weight: < 1.4kg (< 3 lbs.)

### **Output Power**

Up to 10 watts (0-10 watts in 1/10th watt increments on both monopolar and bipolar)

### **Line Voltage**

100-240 VAC

Line Frequency: 50-60 HzMain Current: 320 mA Max.

### **Warranty: FOUR YEARS**

\* Not Pictured



### **ELECTROSURGERY PENCILS**

Most Bovie® electrosurgery pencils utilize a three-pin industry-standard connection and are compatible with most electrosurgical generators. The disposable pencils feature standard pushbutton, rocker switch or foot-control (1-pin connection) activation. The pencils come in a variety of individually packaged sterile kit configurations that include holsters, scratch pads, stainless steel electrodes or Resistick II™ Coated Electrodes. All disposable pencils are packaged individually sterile. Additionally, Bovie® offers two reusable pencils with 3-pin connection for high watt units. The ESPR2 fits standard 3/32″ electrodes and is validated for up to 100 autoclave cycles. The ESPR-EL fits special 1/16 electrodes and is validated for up to 25 autoclave cycles.

		REUS/	ABLE PENCILS		
0 0	ESPR2* Autoclavable reusable pencil, validated to one hundred (100) autoclave cycles, non-sterile - 1/box (electrode not included)			*0 "	
	ESPR-EL* Autoclavable reusable pencil, validated to twenty-five (25) autoclave cycles, 1/16th" collet, non-sterile - 1/box (electrode not included)				g using a standard
		942, DERN	3-button pencil 1 102, DERM 101 - or 1/box		
		DISPOS	SABLE PENCILS		
	ESP1*	ESP1 PUS CONFIGU	H-BUTTON Rations		CK II™ PENCIL Irations
	ESP1T*	ESP1 ESP1H	pencil, blade - 50/box pencil, blade, holster -	ESP1T	disposable pencil, ES01T coated blade - 50/box
	ESP1TN*	ESP1HS ESP1N ESP1HN	40/box pencil, blade, holster, scratch pad - 40/box pencil, needle - 50/box pencil, needle, holster - 40/box	ESP1TN	disposable pencil, ES02T coated blade - 50/box
		ESP6 ROC	KER SWITCH RATIONS		CK II PENCIL Irations
	ESP6	ESP6H	pencil, blade - 50/box pencil, blade, holster - 40/box	ESP6T	disposable pencil, ES01T coated blade - 50/box
		ESP6HS	pencil, blade, holster, scratch pad - 40/box		
	ESP6T*	ESP6N ESP6HN	pencil, needle - 50/box pencil, needle, holster - 40/box		
		ESP7 F00	T-CONTROL CONFIGURATION	ONS	† Pencils using a single pin
	ESP7 <sup>†</sup>	ESP7 ESP7H	pencil, blade - 50/box pencil, blade, holster - 40/box	(	connector

### **ADDITIONAL ACCESSORIES**

Disposable handpiece sheaths are used in conjunction with electrosurgery pencils to provide a cost-effective alternative to continuous sterilization. The use of handpiece sheaths significantly increases the life of the reusable pencils.



<sup>\*</sup> Pencil not included


### **SUCTION COAGULATORS**

Suction coagulators allow fluid evacuation with or without simultaneous coagulation. Suction coagulators are compatible with Bovie® generators, as well as most competitors' generators. The devices feature three different diameters, 6" flexible shaft and hand or foot activation. The single-use devices are individually packaged sterile, 10 per box.



### HANDSWITCH SUCTION COAGULATORS

Catalog Number	Description	Tube Diameter (French Size)	Cable Length
SCH08	Disposable	08 Fr	3m
SCH10	Disposable	10 Fr	3m
SCH12	Disposable	12 Fr	3m

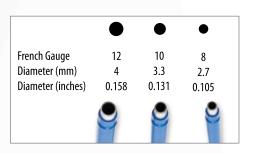




### FOOTSWITCH SUCTION COAGULATORS


Catalog Number	Description	Tube Diameter (French Size)	Cable Length
SCF08	Disposable	08 Fr	3m
SCF10	Disposable	10 Fr	3m
SCF12	Disposable	12 Fr	3m

### **ELECTRODES**

All electrodes utilize the standard 3/32" (2.3mm) stainless steel shaft. The molded electrodes feature safety grip insulators combining patient and user safety with easy insertion into, and removal from, the surgical pencil. Electrodes are manufactured to the highest standards. They come individually packaged sterile or in non-sterile options.\*

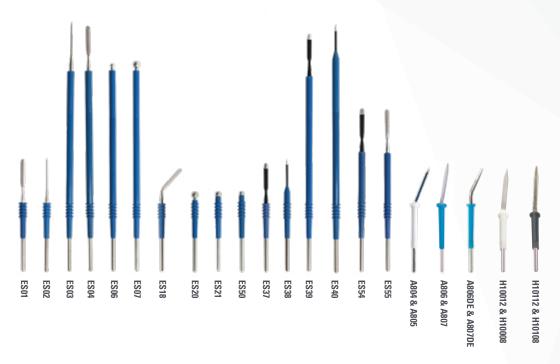



<sup>\*</sup> non-sterile and not packaged individually

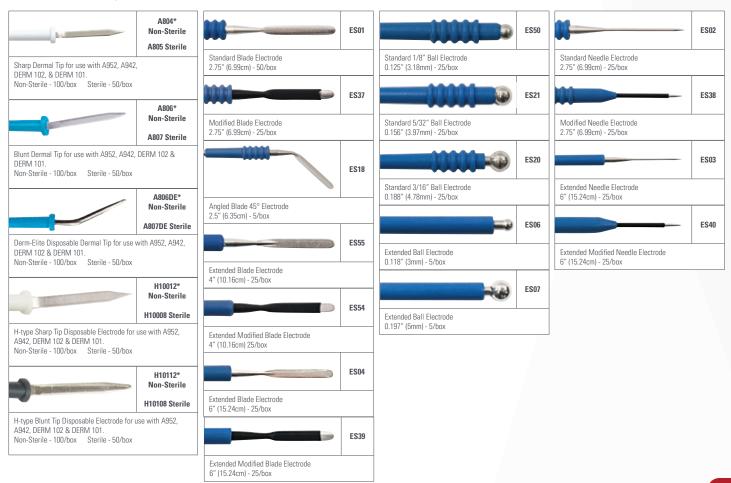


### RESISTICK II™ COATED ELECTRODES

Resistick II<sup>™</sup> coated blade, ball and needle electrodes are expertly coated with PTFE (polytetrafluoroethylene). Resisting eschar build-up, the electrode moves more effectively through tissue, increasing efficiency during surgical procedures. Simply use a damp sponge or cloth to clean any debris off the Resistick II<sup>™</sup> electrodes.

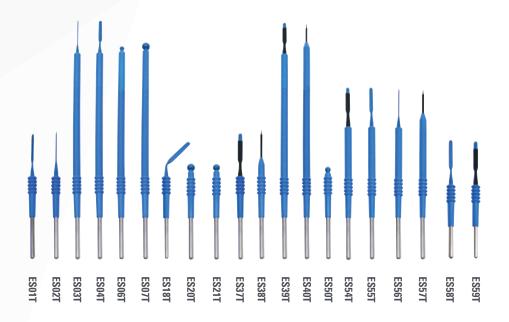

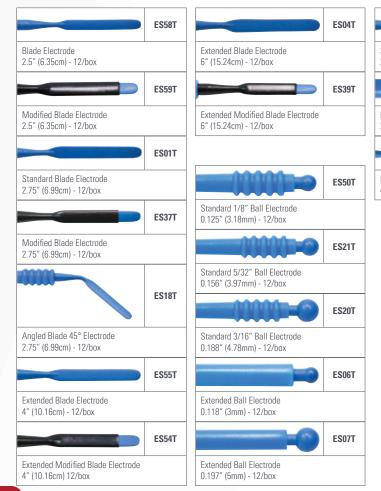


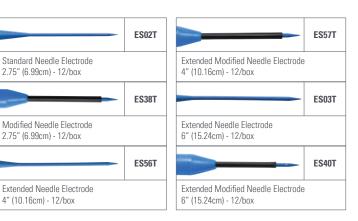

### MICRODISSECTION NEEDLE ELECTRODES

Olsen® and Bovie® Microdissection needles come in a wide array of configurations for a variety of surgical procedures. The ultra sharp design coupled with the finest tungsten wire allows for optimal control and better outcomes for any procedure requiring precise tissue dissection. The fine needle tip concentrates the electrosurgical energy, allowing for low-power settings and less risk of thermal spread.




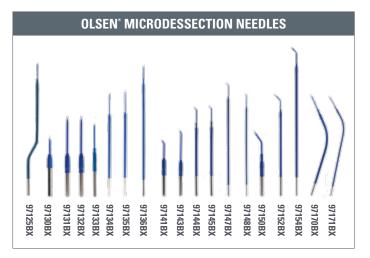


SUPERCUT TUNGSTEN NEEDLES					
	ES62R	ES60R	ES63	ES62	ES61 ES60
		1	1	1	11
			1	1	
	/ '	/			
RE	USABLE SUPER	CUT		SINGL	E-USE SUPERCUT
TL	JNGSTEN NEEDI	.ES		TUN	GSTEN NEEDLES
ES60R	needle super fine 30	cm,	ES		needle super fine 3cm,
ES62R	needle super fine 4.	.5cm,	ES		modified needle super fine 3
CIN	GLE-USE SUPER	CUT			cm, sterile - 5/box
	JNGSTEN NEEDI		F06		U C 4.5
ES61HS	modified fine needle		ES	-	needle super fine 4.5cm, sterile - 5/hox
E901H2	Heatshrink, sterile - !				
	The state of the s		ES		modified needle super fine
-				4	1.5cm, sterile - 5/box




Wire, 14cm (5-1/2 in)

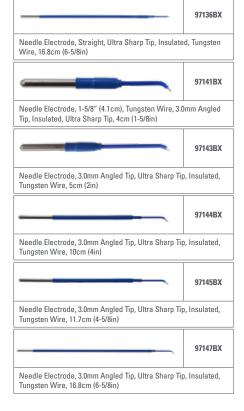

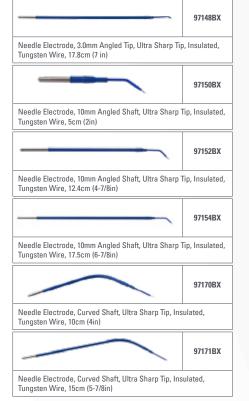

### LAPAROSCOPIC ELECTRODES

Our Bovie\* laparoscopic electrodes are used in minimally invasive procedures performed in outpatient facilities, surgery centers and operating rooms. The electrodes are 33cm long, with a 3/32" shaft, for insertion into standard electrosurgery pencils. Each electrode adaptor is also packaged with an optional 4mm shaft adaptor for connection into 4mm active cables. The laparoscopic electrode tip styles come in a variety of popular configurations with standard stainless steel finish.

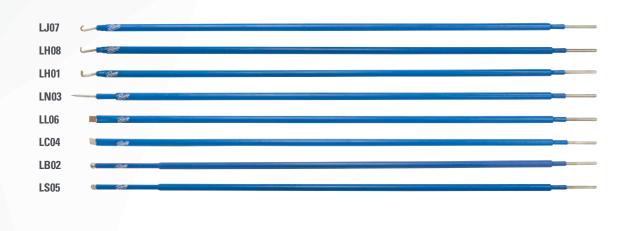


4mm Adaptor






### **DISPOSABLE LOOP ELECTRODES**

All Bovie LLETZ Loop and Square electrodes feature tungsten wire for superior shape and integrity throughout the excision procedure. The electrodes use a standard 3/32" (2.3mm) shaft. All molded loops are packaged in a patented protective tray to prevent damage during shipping.

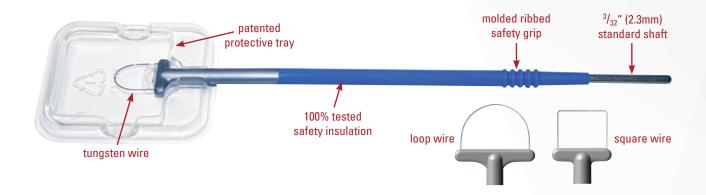

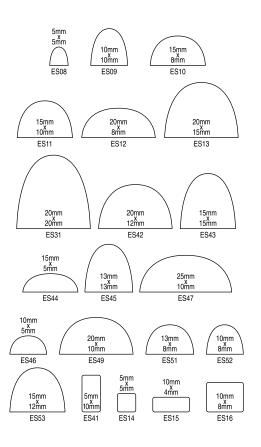

	LOOPS	
ES08 5mm x 5mm loop 5/box	<b>ES31</b> 20mm x 20mm loop 5/box	<b>ES47</b> 25mm x 10mm loop 5/box
<b>ES09</b> 10mm x 10mm loop 5/box	<b>ES42</b> 20mm x 12mm loop 5/box	<b>ES49</b> 20mm x 10mm loop 5/box
<b>ES10</b> 15mm x 8mm loop 5/box	ES43 15 mm x 15mm loop 5/box	<b>ES51</b> 13mm x 8mm loop 5/box
<b>ES11</b> 15mm x 10mm loop 5/box	<b>ES44</b> 15mm x 5mm loop 5/box	<b>ES52</b> 10mm x 8mm loop 5/box
<b>ES12</b> 20mm x 8mm loop 5/box	<b>ES45</b> 13mm x 13mm loop 5/box	<b>ES53</b> 15mm x 12mm loop 5/box
<b>ES13</b> 20mm x 15mm loop 5/box	<b>ES46</b> 10mm x 5mm loop 5/box	

SQUARES				
ES41	ES14	ES15	ES16	
5mm x 10mm square - 5/box	5mm x 5mm square - 5/box	10mm x 4mm square - 5/box	10mm x 8mm square - 5/box	

### **SPECIALTY ELECTRODES**

Symmetry Surgical® manufactures an expansive line of Bovie® specialty electrodes for a variety of applications including Dermal Loop Electrodes for dermatologic procedures, Arthroscopic electrodes and SuperCut™ Needle Electrodes.



### **DISPOSABLE SHORT TUNGSTEN .008" WIRE LOOPS**

**ES22-8** 1/4" (6.35mm) short shaft loop, sterile - 5/box

**ES23-8** 5/16" (7.94mm) short shaft loop, sterile - 5/box

**ES24-8** 3/8" (9.53mm) short shaft loop, sterile - 5/box

**ES25-8** 1/2" (12.7mm) short shaft loop, sterile - 5/box

**ES26-8** 3/4" (19.05mm) short shaft loop, sterile - 5/box

### **DISPOSABLE SHORT TUNGSTEN .015" WIRE LOOPS**

ES22 1/4" (6.35mm) short shaft loop, sterile - 5/box

ES23 5/16" (7.94mm) short shaft loop, sterile - 5/box

ES24 3/8" (9.53mm) short shaft loop, sterile - 5/box

**ES25** 1/2" (12.7mm) short shaft loop, sterile - 5/box

ES26 3/4" (19.05mm) short shaft loop, sterile - 5/box

### **REUSABLE SHORT TUNGSTEN .015" WIRE LOOPS**

**ES22R** 1/4" (6.35mm) short shaft loop, non-sterile - 1/box

ES23R 5/16" (7.94mm) short shaft loop, non-sterile - 1/box

**ES24R** 3/8" (9.53mm) short shaft loop, non-sterile - 1/box

ES25R 1/2" (12.7mm) short shaft loop, non-sterile - 1/box

ES26R 3/4" (19.05mm) short shaft loop, non-sterile - 1/box

### ARTHROSCOPIC ELECTRODES



### **DISPOSABLE ARTHROSCOPIC ELECTRODES**

**AR00** hook 90, sterile - 5/box


AR01 hook 45, sterile - 5/box

AR02 acromioplasty, sterile - 5/box

**AR03** menisectomy, sterile - 5/box

### **REUSABLE ELECTRODES**

All reusable electrodes utilize the standard 3/32" (2.3mm) stainless steel shaft. The molded electrodes feature safety grip insulators combining patient and user safety with easy insertion and removal from the surgical pencil. Reusable electrodes are designed for five -25 autoclave cycles. For care and handling, refer to the Instructions for Use in the package insert.

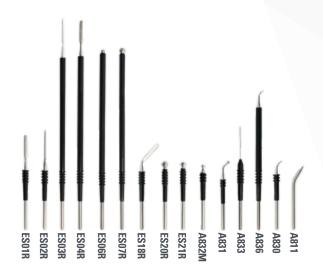









### **REUSABLE LOOPS/SQUARES, ELECTRODES**

All Bovie® LLETZ Loop and Square electrodes feature tungsten wire for superior shape and integrity throughout the excision procedure. The electrodes use a standard 3/32" (2.3mm) shaft. All molded loops are packaged in a patented protective tray to prevent damage during shipping. Electrodes are designed for five autoclave cycles. For care and handling refer to the Instructions for Use in the package insert.

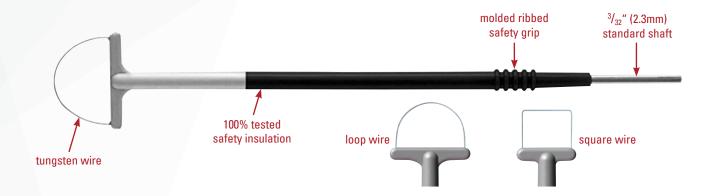

	LOOPS	
ES08R 5mm x 5mm loop 1/box	ES31R 20mm x 20mm loop 1/box	<b>ES47R</b> 25mm x 10mm loop 1/box
<b>ES09R</b> 10mm x 10mm loop 1/box	<b>ES42R</b> 20mm x 12mm loop 1/box	<b>ES49R</b> 20mm x 10mm loop 1/box
ES10R 15mm x 8mm loop 1/box	<b>ES43R</b> 15mm x 15mm loop 1/box	<b>ES51R</b> 13mm x 8mm loop 1/box
<b>ES11R</b> 15mm x 10mm loop 1/box	<b>ES44R</b> 15mm x 5mm loop 1/box	<b>ES52R</b> 10mm x 8mm loop 1/box
ES12R 20mm x 8mm loop 1/box	<b>ES45R</b> 13mm x 13mm loop 1/box	<b>ES53R</b> 15mm x 12mm loop 1/box
ES13R 20mm x 15mm loop 1/box	ES46R 10mm x 5mm loop 1/box	

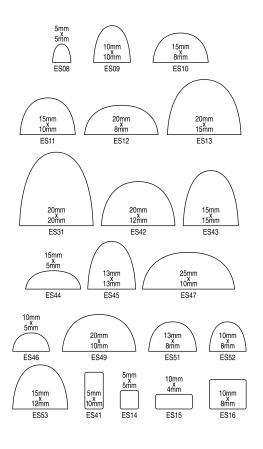

SQUARES			
ES41R	ES14R	ES15R	ES16R
5mm x 10mm square - 1/box	5mm x 5mm square - 1/box	10mm x 4mm square - 1/box	10mm x 8mm square - 1/box

### **REUSABLE BIPOLAR FORCEPS**

Reusable Bipolar Forceps are available in micro (.5mm), 1.0mm, and 2.0mm tip size and in lengths ranging from 3.5" to 8". The forceps are guaranteed for twenty autoclave cycles. For care and handling, refer to the Instructions for Use in the package insert. Reusable bipolar cords, with quality silicone construction are available for all Bovie® generators.

<b>A820</b> 7.25" (184.15mm) bayonet, fine (1.0mm), smooth - 1/box	<b>A821</b> 8" (203.2mm) bayonet, medium (2.0mm), smooth - 1/box	<b>A822</b> 7" (177.8mm) straight, fine (1.0mm), smooth - 1/box
-		
<b>A823</b> 7" (177.8mm) curved, fine (1.0mm), smooth - 1/box	<b>A824</b> 4.5" (114.3mm) curved, micro (.5mm) - 1/box	<b>A825</b> 4.75" (120.65mm) straight, medium (2.0mm), smooth - 1/box
<b>A826</b> 4.5" (114.3mm) straight, micro (.5mm), smooth - 1/box	<b>A840</b> 4" (101.6mm) Tenzel, micro (.5mm), coated - 1/box	A841 3.5" (88.9mm) McPherson, curved, micro (.5mm), coated - 1/box
<b>A842</b> 3.5" (88.9mm) McPherson, straight, micro (.5mm), coated - 1/box	<b>A843</b> 4" (101.6mm) Tenzel, micro (.5mm), uncoated - 1/box	A844 3.5" (88.9mm) McPherson, curved, micro (.5mm), uncoated - 1/box
		➤ Bipolar forceps require an electrosurgery footswitch for
<b>A845</b> 3.5" (88.9mm) McPherson, straight, micro (.5mm), uncoated - 1/box	<b>A827V</b> forceps cord - 1/pouch fits A3350, A2350, A1250S, A952, A942,	activation. See footswitch page for options.

DERM 102, & DERM 101

### **OLSEN® MIDAS TOUCH® BIPOLAR FORCEPS**

The Olsen® Midas Touch® Non-Stick Coating maintains a smooth, highly conductive surface, enabling the use of lower power settings which may reduce tissue damage caused by thermal spread, eschar buildup and the production of toxic smoke.





### **OLSEN® ACCESSORIES**



### **OLSEN® BIPOLAR FORCEPS**

		A Commission of the Commission
<b>20-0421</b> Bayonet, Insulated, Irrigating, Smooth, Straight, 1.0mm Tips, 210mm (8-1/4in) - 10/box	<b>20-0421R</b> Kit, Bayonet, Rotary, Insulated, Irrigating, Smooth, Straight, 1.0mm Tips, 210mm (8-1/4in) - 10/box	<b>20-0431</b> Bayonet, Insulated, Irrigating, Smooth, Straight, 1.5mm Tips, 210mm (8-1/4in)- 10/box
20-1060* Jewelers, Smooth, Straight Micro Tips, 114mm (4-1/2in) - 10/box	20-1060I* Jewelers, Insulated, Smooth, Straight Micro Tips, 114mm (4-1/2in) - 10/box	20-1080* Adson, Serrated, Straight, 1.0mm Tips, 127mm (5in) 10/box
20-1080I* Adson, Insulated, Serrated, Straight, 1.0mm Tips, 127mm (5in) - 10/box	20-1082I* Adson, Insulated, Smooth, Straight, 1.0mm Tips, 127mm (5in) - 10/box	<b>20-1090*</b> Semkin, Smooth, Straight, 0.75mm Tips, 146mm (5-3/4in) - 10/box
<b>20-1190KI</b> Kit, Semkin, Insulated, Smooth, Straight, 1.0mm Tips, 140mm (5-1/2in) - 10/box	<b>20-1360I</b> Bayonet, Insulated, Smooth, Straight, 0.75mm Tips, 181mm (7-1/8in) - 10/box	<b>20-1370*</b> Bayonet, Smooth, Straight, 1.5mm Tips, 181mm (7-1/8in) - 10/box
20-13701* Bayonet, Insulated, Smooth, Straight, 1.5mm Tips, 181mm (7-1/8in) - 10/box	<b>20-1380K</b> Kit, Bayonet, Smooth, Straight, 1.5mm Tips, 191mm (7-1/2in) - 10/box	<b>20-1380I</b> Bayonet, Insulated, Smooth, Straight, 1.5mm Tips, 19mm (7-1/2in) - 10/box
<b>20-1400KI</b> Kit, Bayonet, Insulated, Smooth, Straight, 0.75mm Tips, 210mm (8-1/4in) - 10/box	<b>20-1410</b> Bayonet, Smooth, Straight, 1.5mm Tips, 210mm (8-1/4in) - 10/box	20-1410I* Bayonet, Insulated, Smooth, Straight, 1.5mm Tips, 210mm (8-1/4in) - 10/box
<b>20-1415KI</b> Kit, Bayonet, Insulated, Smooth, Straight, 2.0mm Tips, 210mm (8-1/4in) - 10/box	<b>20-1440I*</b> Bayonet, Insulated, Smooth, Straight, 1.0mm Tips, 210mm (8-1/4in) - 10/box	<b>20-1480I</b> Bayonet, Insulated, Smooth, Straight, 1.0mm Tips, 267mm (10-1/2in) - 10/box

<sup>\*</sup>Available as kit - Add K to PART CODE (Kit contains Forceps, Olsen® Standard Twin Pin Cord and Surgiwipe). Additional Olsen® Bipolar Foreps available. Please contact customer service at 1-800-251-3000.



### **STANDS, FOOTSWITCHES & WALL MOUNT KIT**

When footswitch activation is preferred, Bovie® generators are compatible with the footswitches below. If mobility is a concern, select from a variety of mobile stands.







### RETURN ELECTRODE, ACTIVE CABLE, CORD AND ADAPTOR

Return electrodes are used to safely remove electrosurgical current from the patient. Bovie disposable, single-use return electrodes are manufactured with an adhesive hydrogel for secure patient application as well as reusable return electrode cables and disposable pads to the ESU. Bovie recommends always using a split pad for maximum patient safety. In addition, Bovie provides active cables and an adaptor for the connection of a variety of compatible electrosurgery accessories.

### REUSABLE RETURN ELECTRODE PLATES AND CORDS



**A802EU** metal return electrode plate for A942, DERM 102, & DERM 101 - 1/pouch



A1204P replacement metal return electrode plate for use with A1204C and A1254C - 1/pouch



**A1254C** return electrode cord for A3350, A2350, A1250S, & A952 - 1/pouch



A1252C return electrode cord for ESRE and ESRS to A3350, A2350, A1250S, & A952 - 1/pouch

### **DISPOSABLE RETURN ELECTRODE PADS**



**A1202 or ESRS** solid return electrode A1202 - 5/pack ESRS - 50/box



**ESRSC** solid adult return electrode w/ 2.8m cable - 50/box



**ESRE** split adult return electrode - 50/box



**ESREC** split adult return electrode w/ 2.8m cable - 50/box

FOOT-CONTROL ADAPTOR		PENCIL ADAPTORS
A1205A adaptor for connecting footswitch pencil for the DERM 102, A942, and A952	A1255A adaptor plug for connecting footswitching pencil for the A3350, A2350, A1250S, and most ESUs	A905B pencil adaptor for 1/32 in (0.8 mm) electrode - fits most ESU pencils



### **BOVIE® SMOKE SHARK® II SMOKE EVACUATOR**

The Bovie Smoke Shark II is designed for ease of use and effective smoke and particle removal at the operative site. The Bovie Smoke Shark II is designed for use with all electrosurgical generators and effectively controls smoke plume from lasers and other forms of advanced energy.



### STYLE # SE02

### **▼** Products Included

**SE02** (1) Smoke Shark® II with pneumatic footswitch

**SF35** (1) 35-hour filter

110 VAC Hospital-grade power cord\* User's guide (MC-23088)\*

### **DIMENSIONS**

Height: 15.2cm (6") Width: 27.9cm (11") Depth: 39.4cm (15.5") Weight: 4.3kg (10 lbs.)

### **LINE VOLTAGE**

120 VAC

### LINE FREQUENCY

 $50-60\;\mathrm{Hz}$ 

### **WARRANTY**

Two years

### **TYPICAL SETUPS**

### **ELECTROSURGERY**

**SE02** Smoke Shark unit with pneumatic footswitch (FSSEP)

**SF35** long-life 35-hour filter

**786T** † 7/8" x 6' (22.3mm x 182.9cm ) tube case of 24

**SEAS**\* † arm stand

### **OB/GYN**

**SE02** Smoke Shark unit with pneumatic footswitch (FSSEP)

**SF35** long-life 35-hour filter

**786**T <sup>†</sup> 7/8" x 6' (22.3mm X 182.9cm) tube case of 24 or

**SERF\*** † reducer fitting 7/8" to 1/4" x 24" (22.225mm to 6.35mm x 609.6mm) [fits optional vaginal

speculum]

### **LASER**

**SE02** Smoke Shark unit with pneumatic footswitch (FSSEP)

**SF35** long-life 35-hour filter

**786T** † 7/8" x 6' (22.3mm X 182.9cm ) tube case of 24 or

**SELW\*** † 7/8" x 8"(22.3mm X 182.9cm) laser resistant wand 3/pak


SEAS\* † arm stand



- \* Not Pictured
- <sup>†</sup> Sold Separately

### **BOVIE® SMOKE SHARK® III SMOKE EVACUATOR**

The Bovie® Smoke Shark® III incorporates powerful and effective smoke evacuation to reduce the biological, physical and chemical hazards present in electrosurgical smoke plume. It is designed as a smaller device to uniquely address the limited space and shorter procedures performed in the physician's office.

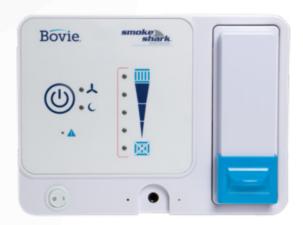

### **STYLE # SE03-120**

### **▼** Products Included

**SE03-120** (1) Smoke Shark\* III **SF06** (1) 6-hour filter **SS95-1** (1) SharkSkin adaptor

SERS2 (1) Smart Shark Remote Activator

Switch

**SE03-M** (1) Wall mount kit 110 VAC Hospital-grade power cord\* User's guide **(P000027200)**\*

### **DIMENSIONS**

Height: 18cm (7.1") Width: 23.9cm (9.4") Depth: 15.5cm (6.1") Weight: 2.8kg (6.25 lbs.)

### **LINE VOLTAGE**

120 VAC

### LINE FREQUENCY

50 - 60 Hz

### **WARRANTY**

Two years

	PRODUCTS	INCLUDED	
SF06 Smoke Shark* III Multi-procedure 6-hour filter	SS95 SharkSkin SharkSkin Electrosurgery Pencil Adaptor with 3/8" x 10' tubing	SERS2 Bovie® Smart Shark Remote Activator Switch	SE03-M Wall mount kit

### **SMOKE EVACUATION PENCILS & ADAPTORS**



### **Smoke Evacuation Pencil**

- Optimal Visibility At Tip
- Powerful Protection From Smoke Plume Hazards
- Removes Up To 55% More Smoke Plume Than Leading Competitors
- Low Profile-Ergonomic Design
- **Economical Solution For Smoke Evacuation**
- Universal Compatibility

**Easy to use CUT and COAG controls** 

**Optimal visibility** 

**Coated Blade Electrode** 

**Ergonomic Integrated Wire** 

Bovie.

**Ergonomic Shape** 

**Tubing Swivels for Ergonomic Movement** 

**Ribbed Body for Improved Grip** 

UP TO STORY WITH ORCA®

**ISEP1000** 

Bovie® ORCA®, Smoke Evacuation Pencil, 10ft Tubing, Sterile, 10/bx



### SS95 SharkSkin

SharkSkin Electrosurgery Pencil Adaptor with 3/8" x 10' tubing, 20/box



Smoke Evacuation Pencil Adaptor and Tubing-3/8" x 4' Hose with 7/8" to 6' Tube, 10 per case



Smoke Evacuation Pencil Adaptor 3/8" x 1' Hose with 3/8" to 3/8" Adaptor, 25 per case

### **SMOKE EVACUATION ACCESSORIES**

Bovie	SF35 long-life 35-hour filter for SE02 model , 1/box		SETW <sup>†</sup> 7/8" x 6' tube with wand & tip, non-sterile, 24/box
	SF18 <sup>†</sup> 18-hour filter for SE01 model, 1/box		SETWS <sup>†</sup> 7/8" x 6' tube with wand & tip, sterile, 24/box
	SF06 Multi-procedure 6-hour filter for SE03 model, 2/box		SEVL <sup>†</sup> 7/8" x 6' tubing reduced to 3/8" for ValleyLab & IEC pencil, sterile, 5/box
Pencil Sold Separately †	SEPA <sup>†</sup> Smoke Evacuation Pencil Adaptor 3/8" x 1' Hose with 3/8" to 3/8" Adaptor, 25 per case		SERS2 <sup>†</sup> Bovie <sup>®</sup> Smart Shark Remote Activator Switch, 1/box
Pencil Sold Separately †	SEPAT † Smoke Evacuation Pencil Adaptor and Tubing- 3/8" x 4' Hose with 7/8" to 6' Tube, 10 per case		SERF <sup>†</sup> reducer fitting 7/8" to 1/4" x 24" 10/box [fits optional vaginal speculum]
	SS95 SharkSkin SharkSkin Electrosurgery Pencil Adaptor with 3/8" x 10' tubing, 20/ box		SERFS † reducer fitting, sterile 10/box [fits optional vaginal speculum]
	FSSEP Pneumatic footswitch included with SE02, 1/box		SELW <sup>†</sup> laser-resistant wand 7/8" x 8", non-sterile, 3/box
	<b>786T</b> <sup>†</sup> 7/8" x 6' tube, non-sterile, 24/box	-1-0	SEAS <sup>†</sup> arm stand, 1/box
	<b>786TS</b> † 7/8" x 6' tube, sterile, 24/box		

<sup>†</sup> Sold Separately

### **SMOKE EVACUATION**

- As the concern for potential dangers of electrosurgically-generated smoke
  plume increases among surgical staff, regulations regarding this matter
  are also on the rise. Some states have passed legislation requiring smoke
  evacuation policies in facilities that use electrosurgery, while others are
  most likely to follow suit.
- Several healthcare societies are recommending the use of smoke evacuation when performing electrosurgical procedures, including the American College of Surgeons (ACS), the Society of American Gastrointestinal and Endoscopic Surgeons (SAGES), and the Association of Perioperative Registered Nurses (AORN).
- An estimated 500,000 healthcare professionals are exposed to surgical smoke proven to carry viruses, carbon monoxide and other chemical and biological hazards¹. Bovie® Smoke Shark™ smoke evacuation systems

- minimize exposure to toxic smoke plume during all forms of advanced energy.
- Bovie® Smoke Shark™ smoke evacuators are compatible with most electrosurgical generators. The units can operate in conjunction with the optional Bovie® Smart Shark Activator (SERS2) for maximized filter life and surgical efficiency.

<sup>1</sup>Barret, Dr. William L., Garber, Dr. Shawn M, "Surgical Smoke - A Review in Literature", Business Briefing: Global Surgery (2004)

### **MEDICAL LIGHTS**

### MI 150 AND MI 300 GOOSENECK EXAM LIGHTING

The best in high-value, space-saving, color-accurate LED exam lights. The MI 300 offers two brightness settings and its mobile base allows it to be easily moved. Both exam lights provide brilliant, color accurate illumination and clarity. All lighting equipment (excluding normal replacement parts) is backed by a three year limited warranty.





Data:	
Light Source	LED
LED Life in Hours	50,000 (average)
Temperature	4,300° Kelvin
Lux (MI 150)	15,000 at 24" (60cm)
	25,000 at 18" (45.7cm)
Lux (MI 300)	30,000 or 15,000 at 24" (60cm)
	50,000 or 25,000 at 18" (45.7cm)
CRI	97
Operating Voltage	100 - 240 VAC 50/60 Hz
Light Hed Input Voltage	24 VDC



Made in the USA.
Only available for sale in the USA.

### MI 150 and MI 300

XLDEM-150 MI 150 Exam Light Connectible Goose Neck w/Fixed Base

(061223)

 ${\tt XLDEM-300} \qquad {\tt MI~300~Exam~Light~Flexible~Goose~Neck~w/Rolling~Base}$ 

(061123)

### **Dimensions:**

Unit Height


MI 150 64.0" (1,625mm) MI 300 66.0" (1,676mm)



### **MI 550 LED EXAM LIGHTING**

The MI 550 light is part of the MI LED Series offered by Bovie\*. The light emits pure white illumination at very low heat with an output rating of 55,000 lux at 4,300 degrees Kelvin color temperature. Engineered with "green" technology, the MI 550 uses significantly less electricity than standard halogen lights. The long-life LEDs are rated at 50,000 hours, virtually eliminating the need for inconvenient bulb changes. All lighting equipment (excluding normal replacement parts) is backed by a three year limited warranty.



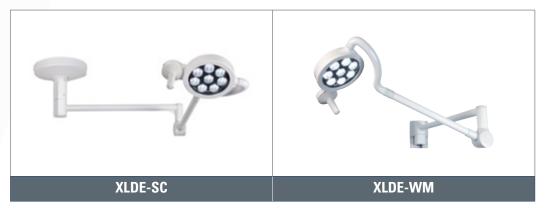




### Dimensions: Head 9.2" (234mm) Dia. x 2.6" (66mm) Deep Arm 54" (1371.6mm) Long Arm (Floor) 21" (533.4mm) Long Floor Stand 74.6" (1895mm) High

Weights:	
Head	6.9lbs (3.1kg)
Arm (Single)	8.3lbs. (21.9kg)
Arm (Dual)	77.1lbs. (35kg)
Wall Bracket Assy.	7.3lbs. (3.3kg)
Floor Stand Assy.	48.4lbs. (21.9kg)
Single Ceiling Mount	145lbs. (65.7kg)
Dual Ceiling Mount	231lbs. (104.7kg)

Data:	
Light Source	LED
LED Life in Hours	50,000 (average)
Color Temperature	4,300° Kelvin
Lux	55,000 Lux at 24" (60cm)
CRI	97
Power Consumption	n 16 Watts
Input Voltage	100 - 240 VAC 50/60 Hz

Rotations:	
Articulating Arm Vertical	+20° Up, -60° Down
Movement	
Articulating Arm Horizontal	≈540°
Movement	
Articulating Arm/Yoke Interface	≈540°
Yoke/Lamp Head Interface	≈300°

MI 550 LE	D Minor Surgery Lights
XLDE-WM	MI 550 Wall Mount: 100-240 VAC, 1.5A, 50/60 Hz, 65W (061322)
XLDE-FM	MI 550 Portable Floor Mode: 100-240 VAC, 1.5A, 50/60 Hz, 65W (061323)
XLDE-SC	MI 550 Single Ceiling: 100-240 VAC, 1.5A, 50/60 Hz, 65W (061324)
XLDE-DC*	MI 550 Dual Ceiling: 100-240 VAC, 3.0A, 50/60 Hz, 130W (061325)
CRST9	9" Ceiling Rod (8' 8" – 10" Mounts) (1001461-9)
CRST20	20" Ceiling Rod (8' 11"— 9' 8" Mounts) (1001461-20)
ECK-30	30" Ceiling Rod (9' 9" – 10' 6" Mounts) (1001461-30)
ECK-40	40" Ceiling Rod (10' 7" – 11' 4" Mounts) (1001461-40)
ECK-50	50" Ceiling Rod (11' 5" – 12' 2" Mounts) (1001461-50)

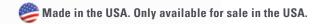



## MI 750 LED EXAM & PROCEDURE LIGHTING

The MI 750 LED Light provides physicians with an economical and energy efficient LED light for examinations and procedures. The MI 750 produces bright, white light for clear and accurate visibility. Designed with unmatched lifetime spring-arm technology, 50,000 hour LEDs, and output rated at 75,000 Lux with a Kelvin color temperature of 4300°, the MI 750 offers a high-quality lighting system that easily aligns with any practice budget. To maintain the sterile field when necessary, the on/off and dimming switches are located on the autoclavable handle. All lighting equipment (excluding normal replacement parts) is backed by a three year limited warranty.







<b>Dimensions</b>	:
Head	14.5" (368mm) Dia. x 3.25" (83mm) Deep
Arm	51" (1295.4mm) Long

18" (457.2mm) Long

74.6" (1895mm) High

Arm (Floor)

Floor Stand

Weights:		
Head	9.0lbs (4.1kg)	
Arm (Single)	48.3lbs. (21.9kg)	
Arm (Dual)	77.1lbs. (35kg)	
Wall Bracket Assy.	7.3lbs. (3.3kg)	
Floor Stand Assy.	48.4lbs. (21.9kg)	
Single Ceiling Mount	145lbs. (65.7kg)	
Dual Ceiling Mount	231lbs. (104.7kg)	

Data:	
Light Source	LED
LED Life in Hours	50,000 (average)
Color Temperature	4,300° Kelvin
Lux	75,000 Lux at 39.4" (1m)
CRI	97
Power Consumption	27 Watts
Input Voltage	100 - 240 VAC 50/60 Hz

Rotations:	
Articulating Arm Vertical	+20° Up, -60° Down
Movement	
Articulating Arm Horizontal	≈540°
Movement	
Articulating Arm/Yoke Interface	≈540°
Yoke/Lamp Head Interface	≈300°

MI 750 LE	D Exam & Procedure Lights
XLDP-WM	MI 750 Wall Mount 100V - 240V (061422)
XLDP-FM	MI 750 Portable Floor Model 100V – 240V (061423)
XLDP-SC	MI 750 Single Ceiling 100V – 240V (061424)
XLDP-DC	MI 750 Dual Ceiling 100V – 240V (061425)
CRST9	9" Ceiling Rod (8' 8"— 10" Mounts) (1001461-9)
CRST20	20" Ceiling Rod (8' 11"— 9' 8" Mounts) (1001461-20)
ECK-30	30" Ceiling Rod (9' 9" – 10' 6" Mounts) (1001461-30)
ECK-40	40" Ceiling Rod (10' 7" – 11' 4" Mounts) (1001461-40)
ECK-50	50" Ceiling Rod (11' 5" – 12' 2" Mounts) (1001461-50)

## MI 1000 LED PROCEDURAL/SURGICAL LIGHTING

The MI 1000 LED light was the first in the MI LED Series offered by Bovie. The light emits pure white illumination at very low heat with an output rating of 100,000 lux at 4,300 degrees Kelvin color temperature. Engineered with "green" technology, the MI 1000 uses significantly less electricity than standard halogen lights. The long-life LEDs are rated at 50,000 hours, reducing the need for inconvenient bulb changes. Five intensity outputs are available and are controlled using the five-stage dimming and on/off switch located on the sterilizable handle. The MI 1000 also has an Endo (Endoscopy) Mode integrated into the lights' dimming capabilities. In this mode, the output produces low level lighting for use during endoscopic surgery. All lighting equipment (excluding normal replacement parts) is backed by a three year limited warranty.







Dimensions	s:
Head	18" (457mm) Dia. x 4" (102mm) Deep
Arm	49" (1244.6mm) Long
Arm (Floor)	16" (406.4mm) Long
Floor Stand	74.6" (1895mm) High

Weight:	
Head	14.6lbs (6.6 kg)
Arm (Single)	48.3lbs. (21.9kg)
Arm (Dual)	77.1lbs. (35kg)
Wall Bracket Assy.	7.3lbs. (3.3kg)
Floor Stand Assy.	48.4lbs. (21.9kg)
Single Ceiling Mount	145lbs. (65.7kg)
Dual Ceiling Mount	231lbs. (104.7kg)

Data:	
Light Source	LED
LED Life in Hours	50,000 (average)
Color Temperature	4,300° Kelvin
Lux	100,000 Lux at 39.4" (1m)
CRI	97
Power Consumption	54 Watts
Input Voltage	100 - 240 VAC 50/60 Hz

						_				
APPL										
Service .	Made	in	the	USA.	Only	available	for	sale	in the	USA.
	iiiuuo		ciio	00/1.	· · · · · · ·	avanabio		Juio		00/1.

Rotations:	
Articulating Arm Vertical	+20° Up, -60° Down
Movement	
Articulating Arm Horizontal	≈540°
Movement	
Articulating Arm/Yoke Interface	≈540°
Yoke/Lamp Head Interface	≈300°

MI 1000 LEI	O Minor Surgery Lights
XLD-WM	MI 1000 Wall Mount 100V - 240V (061522)
XLD-FM	MI 1000 Portable Floor Model 100V - 240V (061523)
XLD-SC	MI 1000 Single Ceiling 100V - 240V (061524)
XLD-DC	MI 1000 Dual Ceiling 100V - 240V (061525)
BV-1000616	MI 1000 Sterilizable Handle (1000616)
CRST9	9" Ceiling Rod (8' 8"— 10" Mounts) (1001461-9)
CRST20	20" Ceiling Rod (8' 11"—9' 8" Mounts) (1001461-20)
ECK-30	30" Ceiling Rod (9' 9"- 10' 6" Mounts) (1001461-30)
ECK-40	40" Ceiling Rod (10' 7" – 11' 4" Mounts) (1001461-40)
ECK-50	50" Ceiling Rod (11' 5" – 12' 2" Mounts) (1001461-50)

## SYSTEM TWO LED SERIES SURGICAL LIGHTING

The System Two LED series presents an unmatched value to surgery centers, specialty practices and outpatient care facilities seeking a surgical lighting system. The energy efficient design provides bright, white, cool and adjustable lighting while consuming half the electricity of its halogen counterpart to offer significant long-term savings for the healthcare facility. The long-life 50,000 hour rated LEDs eliminate procedure room downtime and are output rated at 130,000 Lux with a Kelvin color temperature of 4300°. The System Two LED series is designed with the unique spring-arm technology and has an Endo (Endoscopy) Mode integrated into the lights' dimming capabilities. In this mode, the output produces low level lighting for use during endoscopic surgery. All lighting equipment (excluding normal replacement parts) is backed by a three year limited warranty.



Di	mer	ısio	ns:

Ceiling Mount/Extension Arm Interface Continuous

Articulating/Extension Arm Interface Continuous

Articulating Arm Vertical Movement +15°, -65°

Articulating Arm/Yoke Interface Continuous

Yoke/Lamp Head Interface +/-165° (330° range)

Data:	
Light Source	LED
LED Life in Hours	50,000 (average)
Color Temperature	4,200° Kelvin ±300°K
Lux	130,000 Lux at 39.4" (1m)
CRI	95
Power Consumption	Solo 120W, Duo 240W, Trio 360W
Input Voltage	100 - 240 VAC, 50/60 Hz

Weight:	
Solo Ceiling Mount Assembly	Approximately 73.5lbs. (33.3kg)
Duo Ceiling Mount Assembly	Approximately 104.5lbs. (47.4kg)
Trio Ceiling Mount Assembly	Approximately 155.5lbs (70.5kg)
Light Head Assembly	Approximately 28.1lbs (12.7kg)

Dimensions:	
Ceiling Casting	17.0" (431.8mm) Dia. x 4.7" (119.4mm) Deep
Ceiling Rod	2.5" (63.5mm) Dia. x 9.0"-43.0" (228.6-1092.2mm) Long
Arm Articulated	4.0" (102mm) Dia. (tapered) x 26.0" (660mm) Long
Light Head Assembly	25.5" (648mm) Dia. x 4.75" (120mm) Deep

System Two LED	Surgery Center and OR Lights
XLDS-MA	System Two, Solo, One Monitor Arm (Single Ceiling Mount) (SM1)
XLDS-S2	System Two, Solo, One 130K Lux Light (SD2)
XLDS-S2CR*	System Two, Solo, One In-Light Camera-Ready 130K Lux Light (SD2C)
XLDS-S2MA	System Two, Duo, One 130K Lux Light & One Monitor Arm (SD2M1)
XLDS-S22	System Two, Duo, Two 130K Lux Lights (SD2D2)
XLDS-S22CR*	System Two, Duo, Two 130K Lux Lights (one In-Light Camera-Ready) (SD2CD2)
XLDS-S2MA2	System Two, Trio, One 130K Lux Light & Two Monitor Arms (SD2M1M1)
XLDS-S23MA	System Two, Trio, Two 130K Lux Lights & One Monitor Arm (SD2D2M1)
XLDS-SD2CD2M1*	System Two, Trio, Two 130K Lux Lights (one In-Light Camera-Ready) & One Monitor Arm (SD2CD2M1)
XLDS-S23	System Two, Trio, Three 130K Lux Lights (SD2D2D2)
BV-1000777	In-Light Camera Assembly (1000777)
XLDS-WCK	Wireless wall control kit for System Two lights

<sup>\*</sup>In-Light Camera-Ready System Two lights can be used with or without the camera. To utilize the camera option, items BV-1000777 (In Light Camera Assembly) and XLDS-WCK (Wireless wall control kit for System Two lights) must also be purchased.



#### COLPO-MASTER™ I

Colpo-Master™ colposcopes are a perfect addition to the Bovie® line of gynecologic products. This line of colposcopes offers advanced, 3D, stereoscopic viewing with excellent clarity during gynecologic procedures where millimeters matter. The scopes are available with continuous zoom 4 − 27x magnification. For the EMR ready facility, the Colpo-Master line is available with Trinocular cameraready heads to provide quick and easy visual records of treatments. All Colpo-Master colposcopes are backed by a seven-year limited warranty. Pair it with the Bovie® A1250S-G electrosurgery system with smoke evacuation for a complete GYN LEEP procedure system.



Low center of gravity increases stability, while EZ-Roll locking casters provide smooth gliding movement between exam rooms.

CS-103LED includes 3-legged base as shown. Also available as CS-105LED with 5-legged base.

Colpo-Master™ is a registered trademark of United Products and Instruments, Inc; DBA UNICO

## COLPO-MASTER I ZOOM SUSPENSION ARM COLPOSCOPES

The Colpo-Master is perfect for the GYN performing both exams and procedures, as the swing-arm design allows easy access to the vaginal canal. Colpo-Master I features smooth, continuous zoom 4 – 27x magnification with 20x eyepieces and sharp focus for 3D imaging, and is available with binocular or camera-ready trinocular heads with 298mm working distance. The field of view (FOV) ranges from 58.5mm at 4x to 8.5mm at 27x. All models are available with your choice of 3- or 5-leg sturdy mobile bases. Designed with your viewing comfort in mind, each model has dual diopter adjustments to independently adjust the focus of both left and right eyepieces to eliminate eyestrain. Camera-ready trinocular models allow the user to view through both eyepieces while using the camera.

Colpo-Master I's suspension arm is rock-steady for easy head positioning, and includes locking casters and stainless steel handlebars. It is provided complete with a dust cover, power cord, manual, tool kit and 7-year limited warranty (1-year on the electrical system, excluding bulbs and fuses). Additional accessories are available. Please contact your Bovie® representative for additional information.

Item No.	Description	
CS-103LED	Colpo-Master I Binocular Zoom Suspension-Arm Colposcope	
CS-103T-LED	Colpo-Master I Camera-Ready Zoom Swing-Arm Colposcope	
CS-103T-DV	Swing Arm Digital Video Package, Includes CS-103T- LED and CS-8189 Camera	
CS-103T-HD	Swing Arm HDMI Video Package, Includes CS-103T- LED and CS-8187 HD Camera and Accessories	
CS-103T-HDM	Swing Arm HDMI Video Package, Includes CS-103T- LED and CS-8187 HD Camera and CS-8022 Monitor	
CS-105LED	Same as CS-103LED above, with 5-legged base	
CS-105T-LED	Same as CS-103T-LED above, with 5-legged base	
CS-105T-DV	Swing Arm Digital Video Package, includes CS-105T- LED and CS-8189 Camera	
CS-105T-HD	Swing Arm HDMI Video Package, Includes CS-105T- LED and CS-8187 HD Camera and Accessories	
CS-105T-HDM	Swing Arm HDMI Video Package, Includes CS-105T- LED and CS-8187 HD Camera and CS-8022 Monitor	
CS-8022	21.5" LCD Monitor with HDMI VGA BNC S-Video Inputs	
CS-8187	HD Camera, 1080p with HDMI Output and Accessories	
CS-8189	Digital Video Camera, 14MP, CMOS with USB	

#### **COLPO-MASTER™ II**

#### Colpo-Master™ II

Camera-Ready Trinocular Head with Camera Mount



Colpo-Master  $\!\!\!^{\scriptscriptstyle{\mathrm{T}}}$  is a registered trademark of United Products and Instruments, Inc; DBA UNICO

# COLPO-MASTER II CENTER-POST MOUNT COLPOSCOPES

Colpo-Master II is a compact center-post colposcope with unique, continuous zoom magnification and sharp focusing for crisp 3D imaging.

#### Features

- Full continuous zoom magnification range: 4 27x with 20x eyepieces.
- Binocular or video-ready trinocular heads, with 298mm working distance.
- The field of view (FOV) ranges from 58.5mm at 4x to 8.5mm at 27x.
- 3D Stereoscopic Viewing
- Dual diopters allow user to independently focus each eyepiece to adjust for the natural vision differential between eyes for viewing comfort. Camera-ready Trinocular models allow the user to view through both eyepieces while using the camera.
- Sturdy, 5-legged, easy-to-maneuver base slides in and under exam tables for optimal positioning
- Dual light intensity controls for both white and green LEDs
- Adjustable height center post with smooth suspension system
- HDMI and USB digital camera options available on trinocular version with C-Mount Camera adaptor

It is provided complete with a dust cover, power cord, manual, tool kit and 7-year limited warranty (1-year on the electrical system, excluding bulbs and fuses). Additional accessories are available. Please contact your Bovie® representative for additional information.

Item No.	Description	
CS-205LED	Colpo-Master II Binocular Zoom Center Post Colposcope	
CS-205T-LED	Colpo-Master II Camera-Ready Zoom Center Post Colposcope	
CS-205T-DV	Center Post Digital USB Video Package, Includes CS- 205T-LED and CS-8189 Camera	
CS-205T-HD	Center Post HDMI Video Package, Includes CS-205T- LED and CS-8187 HD Camera and Accessories	
CS-205T-HDM	Center Post HDMI Video Package, Includes CS- 205TLED and CS-8187 HD Camera and CS-8022 Monitor	
CS-8022	21.5" LCD Monitor with HDMI VGA BNC S-Video Inputs	
CS-8187	HD Camera, 1080p with HDMI Output and Accessories	
CS-8189	Digital Video Camera, 14MP, CMOS with USB	

## **OLSEN® VAGINAL SPECULUM WITH SMOKE PORT**



Pictured: 50001

#### OLSEN° VAGINAL SPECULUM WITH SMOKE PORT

The Olsen® Vaginal Speculum is used for performing gynecological pelvic examinations, procedures and Pap Smears. These speculums come in a variety of configurations and include a smoke port for optimal visualization during minor gynecological procedures requiring the use of electrosurgery.


Item No.	Description
50001	Olsen° Speculum, Graves Vaginal Speculum, 3" (7.6cm) X 3/4" (2cm) Blade, Insulated, w/Smoke Tube, 3", 76mm, Reusable
50101	Olsen® Speculum, Graves Vaginal Speculum, 4" (10.2cm) X 1 3/8" (3.5cm) Blade, Insulated, w/Smoke Tube, 4", Reusable
50201	Olsen® Speculum, Graves Vaginal Speculum, 4 1/2" (11.4cm) X 1 1/2" (3.8cm) Blade, Insulated, w/Smoke Tube, 4 1/2", Reusable
50301	Olsen® Speculum, Open Sided Vaginal Speculum, 4 1/2" (11.4cm) X 1 3/8" (3.5cm) Blade, Insulated, w/Smoke Tube, 4 1/2", Reusable



#### BOVIE® ACUDOP™ II DOPPLER SYSTEM

Fetal Dopplers have been used for years to hear and see the fetal heart rate. AcuDop Doppler system provides premium digital sound quality along with breakthrough XD probe technology designed to provide increased sensitivity at shallow depths. Models available in displayed, non-displayed and rechargeable options. Probes 2 MHz, 3 MHz, and waterproof are interchangeable. Convert your OB doppler to a vascular doppler by simply changing the probe to a 5 MHz or a 8 MHz.



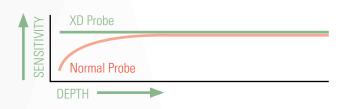

- Unmatched Digital Sound Quality
- Interchangeable OB & Vascular Probes
- Breakthrough XD probes designed to provide increased sensitivity at shallow depths (2-5cm)
- · Cost saving rechargeable models
- Sized to fit comfortably in your hand and pocket
- 5-year Warranty available on ALL Labor and ALL Parts

#### **ACUDOP DOPPLER SYSTEMS**

Style #

AD-770\* AcuDop II 770, Display Unit

**AD-770R** AcuDop II 770R, Display Rechargeable Unit

AD-330\* AcuDop II 330, Non-Display Unit

**AD-330R\*** AcuDop II 330, Non-Display Rechargeable Unit

#### **OB PROBE CONFIGURATIONS**

**OB PROBES:** 

Style #

A2 AcuDop II 2 MHz OB Probe
A3 AcuDop II 3 MHz OB Probe

**WATERPROOF OB PROBES:** 

Style #

A2W AcuDop II 2 MHz Waterproof OB Probe
A3W AcuDop II 3 MHz Waterproof OB Probe



#### **VASCULAR PROBES**


XD technology provides a clear high quality sound to determine patients vascular bloodflow.

#### **VASCULAR PROBE CONFIGURATIONS**

Style #

A5 AcuDop II 5 MHz Vascular Probe
A8 AcuDop II 8 MHz Vascular Probe


\*Not Pictured

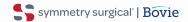




#### **HIGH TEMPERATURE CAUTERIES**


Bovie has been the number one manufacturer of battery-operated cauteries for decades. The safety-enhanced, snap-design cautery is the only battery-operated cautery in the world to address known safety hazards of improper disposal. The improved patent-pending cautery offers a safety-enhanced design, with a recessed button to reduce the potential for inadvertent activation. The cautery units are also modified to enable practitioners to "snap off" the end of the cautery for safe and effective removal of the battery power source.

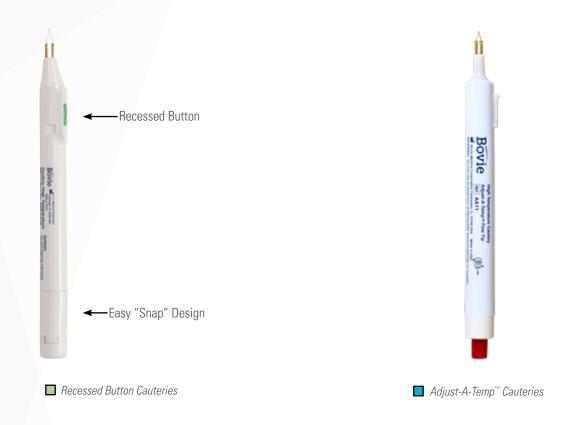


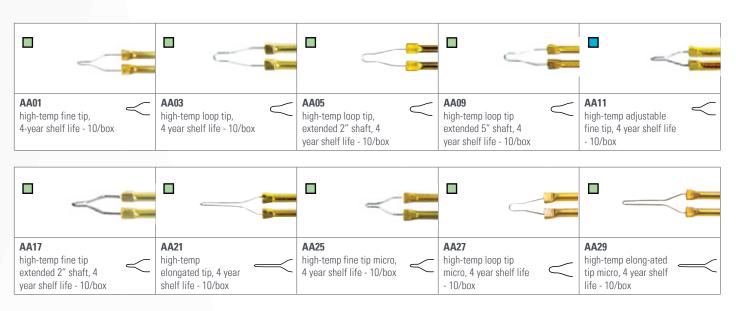

## **LOW-TEMPERATURE CAUTERIES**

Bovie has been the number one manufacturer of battery operated cauteries for decades. The safety-enhanced, snap-design cautery has been enhanced based on valuable customer feedback. The improved cautery offers a new safety-enhanced design with a recessed button to reduce the potential for inadvertent activation. Bovie also modified the cautery units to allow the practitioner to "snap off" the end of the cautery for safe and effective removal of the battery power source. 10 per box.

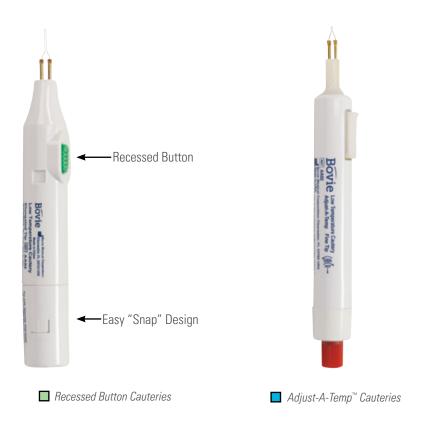

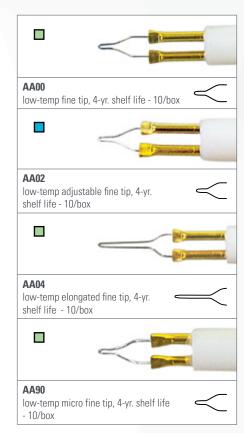



## **CHANGE-A-TIP™ CAUTERIES**

These replaceable battery and replaceable tip cauteries are the next advancement in battery operated cauteries where reusability is a must. These products are designed for the office, ER, and clinical procedures. They can be purchased as the DEL0, DEL1 or DEL2 pre-packaged Change-A-Tip™ cautery and tip kits or individually in high or low temperature models with a variety of tips designed to fit your procedural needs.







Change-A-Tip high-temp handle (includes one H101 non-sterile tip



# Change-A-Tip deluxe **low-temp** cautery kit (includes one low-temperature handle, two H100 non-sterile tips, two H104 non-sterile tips, one AA battery and a foam-lined case)

Recessed Button Cautery

#### DEL1

Change-A-Tip deluxe **high-temp** cautery kit (includes one high-temperature handle, two H101 non-sterile tips, two H103 non-sterile tips, two AA alkaline batteries and a foamlined case) - 1/box



■ Recessed Button Cauteries

#### DEL2

and two AA alkaline batteries) - 1/box

Change-A-Tip deluxe HI-LO cautery kit (includes **one low-temp** handle, **one high-temp** handle, one sterile H100 tip, one sterile H101 tip, one sterile H104 tip, six AA alkaline batteries and a foamlined case) - 1/box

**DELO** 


1/box

## **CAUTERY REPLACEMENT TIPS**

Sterile cautery tips include a sterile drape to cover the non-sterile handle. Tips are individually packaged, 10 per box. The replacement tips are used in conjunction with the Change-A-Tip™ cautery handles as well as the deluxe kits.



Protective Tray

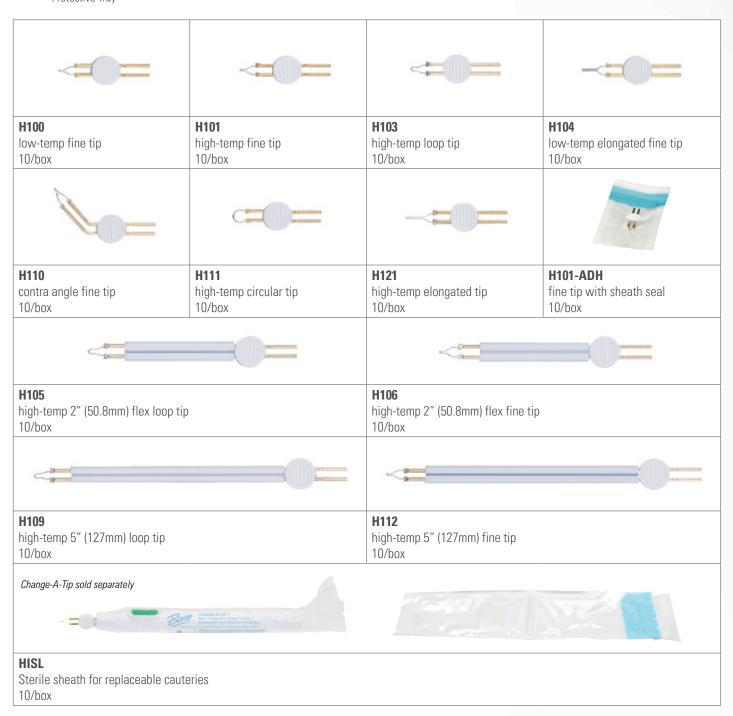



## **FLEXIBLE LIGHTS AND STYLETS**

Bovie\* manufactures two groups of flexible lighting products. The first is available in three different sizes: 5" (127mm), 10" (254mm), and 15" (381mm) (refers to the flexible shaft portion of the light). These lights are for use in procedures where additional light is needed. The second group of flexible lights is our lighted stylets, for difficult intubations.



## **UV WOODS LIGHT, SPECIALTY LIGHTS, PENLIGHTS**



Specialty lights are used in conjunction with fluorescein dye to find foreign matter in the eye.



disorders and rape detection procedures - 1/box

# **PENLIGHT PACKS** Bovie Bovie Diagnostic Medical Penlights High Intensity Medical Penlights · Disposable · Disposable 6303 6666 cobalt penlight - 3/pack doctor's disposable penlight - 6/pack Bovie RN-Pack™ Nurses Examination Light High Intensity Medical Penlights **66RN** 7156 penlight with pupil gauge - 6/pack cobalt filter only - 1/pack



## **DERMASCOUT™ II**

The DermaScout<sup>™</sup> II is an efficient and dependable instrument for visualization of dermal lesions. The dermatoscope is used to view, measure and monitor skin lesions, fungal infections, tumors, warts and other conditions that may affect the patient's epidermis. The hand-held magnifier has 6x magnification, LED illumination, and a dual scale (Eng./Metric) measuring reticle, as well as an independent diopter to allow the user to focus on differing focal points.





**LED** Illumination



**Storage Case** 

#### Catalog #

**BVDS-6-02** – DermaScout™ II Dermatoscope **DS-6-01** – Replacement Battery, Type CR2, 3


#### **Features**

- Magnification with 40mm top lens
- Diopter focus control for sharp focusing
- 40mm field of view (flat field)
- 31mm working distance
- Illuminated with an array of 9 LEDs
- More than 150 hours of battery power
- · Reticle with dual measuring scales: millimeters and inches
- 20mm: 10mm left of center/10mm right of center
- Measure in 0.1mm increments
- Measure in 0.05" increments
- Durable locking transport and storage case
- Battery included (type: CR2, 3V)

#### **OPHTHALMIC PRODUCTS**


Bovie offers a diverse line of ophthalmic products, including an ophthalmic burr, and a clear protective eye bubble for ophthalmic and neurological procedures.



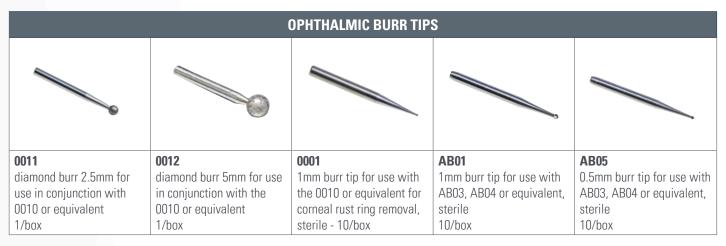

### **NERVE LOCATORS**

Neuro-Pulse<sup>™</sup> nerve locators are designed for identifying exposed motor nerves and reducing the possibility of accidental nerve damage or severance. These battery-powered units are portable and are ideal for head, neck, hand, plastic and facial applications.

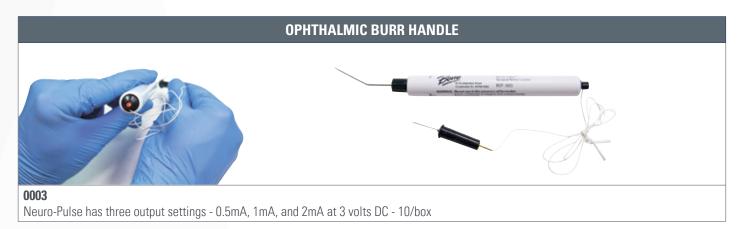

# NOTES: