



With a rich legacy and expansive product portfolio, Symmetry Surgical provides one of the most comprehensive surgical instrument offerings in the marketplace today. We offer multiple instrument styles to meet your performance or procedural needs as well as budget.

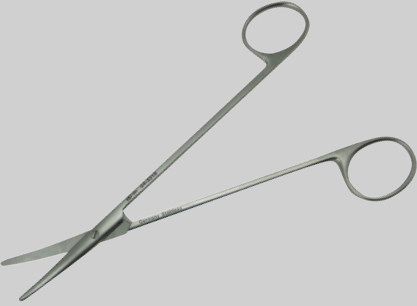
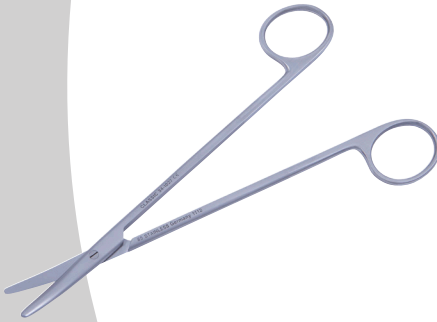
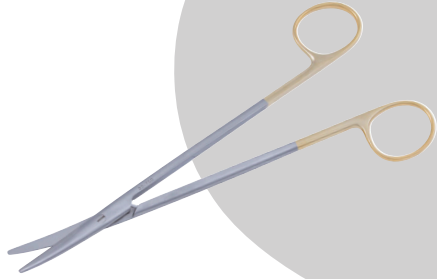



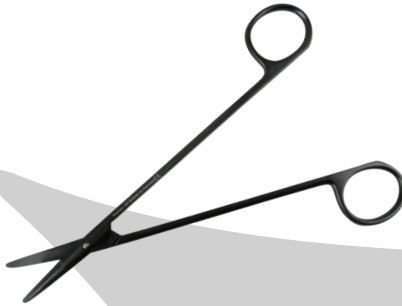
A leader in the surgical instrument market, Symmetry Surgical is 100% focused on surgical instrumentation. Our breath of line allows customers to depend on us for virtually all of their surgical instrument needs and to meet the challenges you face in the OR.

Symmetry Surgical's unparalleled scissors offering represents the strength and diversity of our high-quality instrument portfolio. Whether you need standard stainless steel workhorse scissors or the finest ceramic coated scissors available for the most demanding procedures, Symmetry Surgical has a full spectrum of scissors solutions.

- We offer a wide choice of stainless steel, tungsten carbide, Stellite[®], and ceramic coated scissors.
- Some patterns include serrated lower blades for better cutting action.
- Our Classic[®] and Classic Plus[®] branded scissors offer five and ten year sharpness warranties.

On the reverse of this page, we have listed some options and features for your curved 7 inch Metzenbaum scissors. Learn more about our scissors portfolio as well as other surgical instrumentation at [symmetrysurgical.com](https://www.symmetrysurgical.com).



			
<p>Standard Stainless Steel 04-3318</p>	<p>Classic® Stainless Steel 34-1027</p>	<p>Ultra® Series Tungsten Carbide 32-750 Standard Tips 32-751 Delicate Tips</p>	<p>Classic Plus® Stellite® 36-5016 Standard Tips 36-5033 Delicate Tips</p>
<ul style="list-style-type: none"> • Basic pattern • 30° bevel • Our most economical 7 inch Metzenbaum scissors • German made 	<ul style="list-style-type: none"> • Five year sharpness warranty • Flat screw provides a smooth, clean profile to prevent tissue abrasion • Finish is a unique combination of bright and matte to help distinguish from other instruments • German made 	<ul style="list-style-type: none"> • Durable tungsten carbide cutting edges stay sharp longer than standard stainless steel • Gold ring handles • German made 	<ul style="list-style-type: none"> • 10 year sharpness warranty • Gold ring handles in this pattern designate Stellite® cutting edge inserts • Stellite® is a durable stainless alloy • German made
			
<p>Razor Edge 80-3318</p>	<p>Razor Edge Serrated 86-3318</p>	<p>Ceramic Coated Tungsten Carbide 224.33.68</p>	
<ul style="list-style-type: none"> • Combines the sharpness of a razor-like 10° bevel top and a 30° bevel lower blade • Fine minimally traumatizing cut • The reduced thickness on backside of the cutting blade allows better visibility • Black ring handles • German made 	<ul style="list-style-type: none"> • 10° bevel top blade and serrated lower blade that holds tissue in jaws • Very sharp blades provides minimally traumatic dissection • Sharp enough to dissect by sliding the open blade through tissue • Black ring handles • German made 	<ul style="list-style-type: none"> • Ceramic tungsten carbide cutting edges with Razor Edge on one blade and micro serrations on the other blade • Low glare black finish and gold screw • High-end solution for the most discerning surgeons • Combines the benefits of a Razor Edge serrated style of scissors with the durability of a ceramic coated tungsten carbide edges • German made 	

*Stellite® is a registered trademark of Kennametal Inc.