

Personal Protective Equipment

Isolation & Specialty Gowns

 **aspen**surgical®




Aspen Surgical strives for unmatched quality and unparalleled dependability with brands like Bard-Parker®, Colby™, Precept®, Protek®, Stork®, SurgiSafe®, WriteSite®, and more. Aspen Surgical is a leading provider of surgical essentials, PPE, ultrasound imaging accessories, and labor and delivery disposables. We got here by providing solutions that are as effective as they are trusted. With decades of expertise in the operating rooms, we are not always recognized by name, but our reputation is always at hand. And although we have been supporting medical teams for many years, we never stop evolving, looking ahead, and expanding our reach.

Learn more about what Aspen Surgical has to offer in the pages of this catalog, or online at www.aspensurgical.com



Table of Contents

Isolation Gowns... 6-8	
AAMI Level 2 Full Coverage Gowns pg... 6-8	
Protective Gowns... 9	
Yellow Gowns, Open Back pg... 9	
Blue Film Gowns, Open Back pg... 9	
Specialty Gown... 10	
Chemotherapy Gown pg... 10	

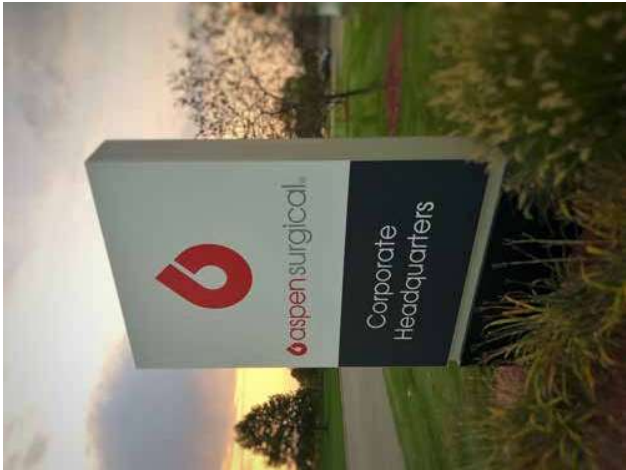

Aspen Surgical Products, Inc. Headquarters, Caledonia, Michigan

Precept® is a registered trademark of Aspen Surgical Products, Inc.
© 2021 Aspen Surgical Products, Inc. All rights reserved.

51177
51198



- Isolation Gown Full Coverage**
Product Details
- Multi-layer SMS material[®] delivers AAMI Level 2 protection and comfort
 - Available in various wrist styles to ensure proper glove donning
 - Full coverage with overlapping back panel and secure neck ties
 - Sonic-welded seams
 - Made in North America
 - Not made with natural rubber latex



Barrier Protection AAMI PB70:2012

Test Method	Test Result	Anticipated Risk of Exposure
AATCC 42 Spray Impact	≤ 1.0g	Low
AATCC 127 Hydrostatic Pressure	≥ 20cm	Low

Part No.	Size	Style	Color	Material	Packaging	Qty
51177	Large	Secure neck tie with elastic wrist	Yellow	SMS	10/Bag	100/Case
51198	Universal	Secure neck tie with knit cuffs	Yellow	SMS	10/Bag	100/Case

51192FC
51623FC



- Isolation Gown Over-the-Head Thumb Loop Full Coverage**
Product Details
- Multi-layer SMS or polycoated material[®] delivers AAMI Level 2 protection and comfort
 - Available in various wrist styles to ensure proper glove donning
 - Full coverage with overlapping back panels
 - Secure neck ties
 - Sonic-welded seams
 - Made in North America
 - Not made with natural rubber latex



Barrier Protection AAMI PB70:2012

Test Method	Test Result	Anticipated Risk of Exposure
AATCC 42 Spray Impact	≤ 1.0g	Low
AATCC 127 Hydrostatic Pressure	≥ 20cm	Low

Part No.	Size	Style	Color	Material	Packaging	Qty
51192FC	Universal	OTH thumb loop with tapered wrist	Yellow	SMS	10/Bag	100/Case
51623FC	Universal	OTH thumb loop with elastic wrist	White	Polycoated	10/Bag	100/Case

59979



Isolation Gown Full Coverage

Product Details

- Multi-layer SMS material[®] delivers AAMI Level 2 protection and comfort
- Elastic wrist to ensure proper glove donning
- Full coverage with overlapping back panels
- Tape tab neck closure
- Sonically sealed seams on sleeves, armholes, shoulders, and side flaps
- Made in North America
- Not made with natural rubber latex



Barrier Protection AAMI PB70:2012

Test Method	Test Result	Anticipated Risk of Exposure
AATCC 42 Spray Impact	≤ 1.0g	Low
AATCC 127 Hydrostatic Pressure	≥ 20cm	Low

Part No.	Size	Style	Color	Material	Packaging	Qty
59979	Universal	Tape tab neck closure with elastic wrist	Yellow	SMS	10/Bag	100/Case

51192



Protective Gown Over-the-Head Thumb Loop with Open Back

- Medium weight multi-layer SMS material[®] passes AAMI Level 2 protection and provides comfort
- Over-the-head donning improves compliance
- Thumb loops and elastic wrist ensure proper glove donning
- Sonic-welded seams
- Made in North America
- Not made with natural rubber latex

Part No.	Size	Style	Color	Material	Packaging	Qty
51192	Large	OTH thumb loop with elastic wrist	Yellow	SMS	10/Bag	100/Case

51575



Blue Film Gown Over-the-Head Thumb Loop with Open Back

- Blue polyethylene film material[®] passes ASTM F1671 viral testing providing enhanced barrier protection
- Generous sizing improves comfort
- Over-the-head donning improves compliance
- Tapered wrists and thumb loops ensure proper glove donning
- Thermal impulse seams
- Made in North America
- Not made with natural rubber latex

Part No.	Size	Style	Color	Material	Packaging	Qty
51575	Universal	OTH open back tapered thumb loop	Blue	Polyethylene Film	10/Bag	100/Case

50530



Chemotherapy Gown, Full Coverage

Product Details

- Tested for commonly used chemotherapy drugs per ASTM F7320®
- Polycocated material delivers AAMI Level 2 protection and comfort
- Knit cuffs provide comfort
- Secure neck ties
- Full coverage with overlapping back panels
- Intended to be used in environments where there is a potential hazard in the handling of chemotherapy drugs
- Universal, one size fits most
- Made in North America
- Not made with natural rubber latex



Barrier Protection AAMI PB70:2012

Test Method	Test Result	Anticipated Risk of Exposure
AATCC 42 Spray Impact	≤ 1.0g	Low
AATCC 127 Hydrostatic Pressure	> 20cm	Low

Part No.	Size	Style	Color	Material	Packaging	Qty
50530	Universal	Secure neck tie with knit cuffs	White	Polycoated	10/Bag	100/Case

Notes:

①AATCC 42 – Water Resistance: Spray Impact Penetration Test / AATCC 127 - Water Resistance: Hydrostatic Pressure Test For more information on these test methods, please visit aatcc.org and aami.org.

②ASTM F1671: Standard test method for resistance of material used in protective clothing to penetration by blood-borne pathogens using Phi-174 bacteriophage penetration as a test system. For more information on this test method, please visit astm.org.

Chemotherapy Commonly Tested Drugs

③The material and seams used in these gowns were tested against the following chemotherapy drugs per ASTM F739-20, and no breakthrough was observed during the test duration of 480 minutes based on average standard breakthrough time.

Carmustine (BCNU)	Cisplatin
Cyclophosphamide (Cytoxan)	Dacarbazine
Doxorubicin HCl	Etoposide
Fluorouracil	Methotrexate
Mitomycin C	Paclitaxel
ThioTepa	Vincristine Sulfate

Aspen Surgical makes no representation or warranty that its products will provide absolute protection from contact or exposure to contaminants or blood, or eliminate the risk of contracting diseases of any nature. OSHA requires the employer to evaluate the task and type of exposure anticipated and select the appropriate personal protective equipment (PPE) that reduces the employee's exposure to potential hazards. Aspen Surgical reserves the right to make changes without notice in design, specifications and models.

Precept® is a registered trademark of Aspen Surgical Products, Inc.

©2021 Aspen Surgical Products, Inc. ALL RIGHTS RESERVED.
M-SM-CAT-1062404 rev1 19-JUL-2021 ENG – US

For further information about this product or a service, please contact your local Aspen Surgical representative or visit our webpage:

USA 616.698.7100
Canada 888.364.7004
www.aspensurgical.com

