

EQUASHIELD®

# Vial Adaptor



# Vial Adaptor

## Product Overview

When handling hazardous drugs, safety is paramount. That's why EQUASHIELD® Vial Adaptors snap on securely and are designed with elastomeric membranes that are airtight and leak-free. EQUASHIELD® Vial Adaptors are available in various sizes, accommodating all standard vial neck sizes.

## Product Features

- / Vial Adaptors are designed to firmly snap onto any standard vial, ensuring a secure and precise fit.
- / Vial Adaptors provide an audible click as a form of feedback, letting you know when they're successfully engaged.
- / Equipped with locking telescopic arms, our Vial Adaptors prevent accidental removal once engaged. This feature aids in minimizing the potential for exposure.
- / The self-centering ring on EQUASHIELD® Vial Adaptors allows for centered spiking each time on neck vials, eliminating angle spiking.
- / Elastomeric membrane wipes the internal needle clean after each use, effectively preventing hazardous drug exposure.

## Fluid Pathway

- Polypropylene
- Silicone
- Polyisoprene
- Hydrophobic Barrier

## Online Resources

[Product Page](#)[How to videos](#)[EQ Academy](#)


## Use with

EQUASHIELD® Syringe Unit,  
available in 13,17,20,28,32 mm

## Ordering Information

### Vial Adaptors

**VA-13/2**  
13mm

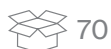

**VA-13C/2**  
13mm Concave




**VA-17/2**  
17mm

Not available  
In the US




**VA-20/2**  
20mm

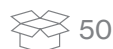

**VA-20C/2**  
20mm Concave



**VA-28/2**  
28mm

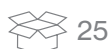

**VA-32C/2**  
32mm

Not available  
In the US



\*Only available outside of the US

[Find a US Sales Representative](#)

[Find a local Distributor](#)

[Contact EQUASHIELD](#)