

EQUASHIELD®

Syringe Units



Syringe Unit

Product Overview

EQUASHIELD®'s Syringe Units, are your frontline defense against contamination from hazardous drug handling. With it's built-in closed pressure equalization system and encapsulated plunger rod, sterility is maintained.

Syringe Units come with a connector that is permanently bonded to the body of the Syringe Unit, significantly reducing the risk of accidental disconnections. Additionally, the needles are shielded, offering an extra layer of protection against accidental needle sticks, further reinforcing safety.

With availability in various sizes and **FDA clearance for full volume use**, Syringe Units are an efficient, easy-to-use solution for your drug handling needs.

Product Features

/ Syringe Units are uniquely encapsulated with a built-in pressure equalization system, preventing plunger rod contamination.

/ The market's only closed-back Syringe Unit, preventing the release of hazardous vapors from the back of the syringe, ensuring superior safety compared to traditional open-back syringes.

/ The welded connector on our Syringe Units drastically reduces the risk of accidental disconnections and decreases the steps to compound a dose.

/ Elastomeric membrane to facilitate dry disconnections, effectively preventing contamination from drug residue.

/ The welded connector on our Syringe Units drastically reduces the risk of accidental disconnections and decreases the steps to compound a dose.

Fluid Pathway

- Polypropylene
- Stainless Steel
- Silicone
- Polyisoprene

Technical Specification

- **Priming Volume:** 0.07 ml for syringe 1 ml to 10 ml and 0.08 ml for all other syringes
- **Needle Size:** 17G (1.5 mm)
- **Needle Length:** 1.2" (30 mm)
Non-Coring pencil point tip: needle never penetrates vial rubber stopper

Online Resources

[Product Page](#)[How to videos](#)[EQ Academy](#)

Use With

EQUASHIELD® Vial Adaptor, Spike Adaptor, Luer Lock Adaptors, Secondary Tubing Set, Luer Lock Y-Site, Luer Lock Adaptor C, Luer lock Adaptor DC, Protective Plug
Available in 1,3,5,10,20,35,60 ml

Ordering Information

Syringe Units

SU-1/2

Syringe Unit
1mL



SU-3/2

Syringe Unit
3mL




SU-5/2

Syringe Unit
5mL



SU-10/2

Syringe Unit
10mL

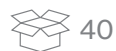


SU-20/2

Syringe Unit
20mL




SU-35/2

Syringe Unit
35mL

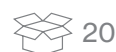


SU-60/2

Syringe Unit
60mL

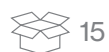

SU-EZ60/2

Syringe Unit
60mL

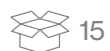

[Find a US Sales Representative](#)

[Find a local Distributor](#)

[Contact EQUASHIELD](#)